

Non-Small Cell Lung Cancer Spheroids to Investigate Factors in Genome Stability



DEPARTMENT OF
Biomedical Engineering
UNIVERSITY OF WISCONSIN-MADISON

Final Report

BME 402: Biomedical Engineering Capstone Design

May 5, 2025

Althys Cao, Emily Rhine, Ana Martinez, Julia Salita, Jayson O'Halloran, Ms. Carley Schwartz, Dr. Gaelen Hess (Department of Biomolecular Chemistry and Center for Human Genomics and Precision Medicine), and Dr. Paul Campagnola (Department of Biomedical Engineering, Peter Tong Department Chair)

Abstract

Lung cancer is the leading cause of cancer death in the US by a large margin, accounting for about 1 in 5 of all cancer deaths according to the American Cancer Society. Approximately 80% of lung cancer cases are classified as non-small cell lung cancer (NSCLC), a condition marked by significant clonal and morphological diversity within a tumor according to Rozenberg et al.. When it comes to *in vitro* models that can accurately mimic the microenvironment of NSCLC, spheroids are an ideal approach. Current two-dimensional models have limited biological relevance to the tumor microenvironment, often failing to replicate the complex cell-cell interactions found *in vivo*. Spheroids can be cultured from a variety of sources including primary isolated tumor cells or continuous cancer cell lines. NSCLC spheroid models can provide a more accurate analysis of intercellular interactions, oncogenic signaling, and anti-cancer drug effects, making them a helpful model in improving the therapies associated with NSCLC. Although our findings from SOX2 real-time quantitative PCR lack the statistical significance to support the claim that there are different genetic expression profiles between 2D and 3D cancer models, a need for an optimized spheroid formation protocol that can be accurately stained for DNA double-stranded breaks ahead of a genome wide CRISPRi screen to find tumor-specific factors that regulate genome stability was established.

Introduction

Non-small cell lung cancer (NSCLC) has a 28% 5-year survival rate, and continues to be one of the most challenging cancers to treat [1]. A major obstacle in treating NSCLC is the continuous evolution of cancer genomes, which leads to frequent mutations, clonal diversity, and resistance to therapeutic treatment. Additionally, factors such as immune evasion, an abundance of extracellular matrix (ECM), and hypoxic conditions further complicate treatment options and are linked to poor patient outcomes [2].

A549 Cell Line

Accurately replicating these characteristics in tumor models is necessary for refining existing therapies and advancing new treatment approaches. For this reason, the A549 cell line was selected to use in a 3D spheroid model. They exhibit adherent properties and key features of lung adenocarcinoma, including epithelial characteristics and high proliferative capacity [2]. Notably, A549 cells mimic the hypoxic core of *in vivo* tumors, which has been shown in NSCLC to drive changes in gene expression, increase oxidative stress, and contribute to chemoresistance, conditions that are poorly represented in traditional 2D culture [2]. The hypoxic core promotes progression by stabilizing HIF-1 α , which drives the expression of genes involved in angiogenesis, metabolic reprogramming, and survival under low-oxygen conditions. Additionally, A549 spheroids exhibit marked clonal heterogeneity, increased cell-cell interactions, and more complex signaling dynamics that mirror the diversity of tumor subpopulations observed in patients [2]. According to Choe et.al, in a SOX2 stemness study, spheroids also showed elevated expression of SOX2, which is a pluripotency marker associated with stem-like traits, enhanced tumor initiation capacity, and resistance to treatment. Taken

together, these characteristics make A549 spheroids a better model at capturing the complexity of NSCLC tumors in vitro compared to 2D monolayer cultures.

SOX2

SOX2 is a transcription factor associated with stemness and tumor progression in NSCLC. It regulates self-renewal, differentiation, and pluripotency, making it a key driver of cancer stem cell (CSC) properties. SOX2 is frequently co-expressed with surface markers such as CD133 and CD44, which are commonly used to identify CSCs. CD133, a marker of tumorigenicity, is often upregulated in SOX2-expressing cells and is linked to enhanced self-renewal and tumor initiation [9]. Similarly, CD44, which plays a role in cell adhesion and migration, is associated with SOX2 in promoting metastasis and resistance to therapies [9]. These surface markers, in conjunction with SOX2, delineate a subpopulation of CSCs that contribute to lung cancer progression, recurrence, and therapeutic resistance, making them a target for gene therapies in NSCLC. Typically SOX2 is more upregulated in spheroid models compared to 2D culture due to microenvironmental cues that mimic in vivo tumor conditions including hypoxia, cell-cell interactions, and mechanical confinement [3].

γ H2AX

The use of γ H2AX staining has been important in understanding the mechanisms of DNA double-strand breaks (DSBs) in 3D spheroid models. γ H2AX is a phosphorylated variant of the histone H2AX, which serves as a sensitive marker for DNA DSBs [4]. The phosphorylation of H2AX occurs rapidly in response to DNA damage, making it an essential marker for monitoring the integrity of the genome [5]. To induce DSBs, etoposide can be introduced into the spheroids. Etoposide is a chemotherapeutic drug that works as a topoisomerase II inhibitor, preventing the re-ligation of DNA strands during replication [6]. It has been seen that compared with 2D cancer

cell cultures, spheroids exhibit a more typical in vivo response to cytotoxic drugs like etoposide [2]. Measuring DNA damage using a 3D γ H2AX staining protocol will allow for better understanding of how disruptions in DSBs repair contribute to NSCLC progression, providing a greater understanding of how 2D and 3D environments vary in cancer modeling.

Materials and Methods

Cell Line and Maintenance

WT A549 cell line was obtained frozen from Hess Lab, which was originally purchased from American Type Culture Collection (ATCC). Cells were thawed by slow mixing of complete DMEM media, which was DMEM media (4.5g/L D-Glucose, L-Glutamine Gibco™ DMEM, Thermo Fisher Scientific, USA, 11965118) supplemented with 1% v/v Gibco™ Penicillin-Streptomycin (Thermo Fisher Scientific, USA, 15070063) and 10% v/v Gibco™ Fetal Bovine Serum (FBS) (Thermo Fisher Scientific, USA, A5256701). Once thawed, cells were pelleted and resuspended in complete DMEM.

Cells were cultured and passaged on T75 cell culture flasks (VWR, USA, 10062-860). Cell cultures were incubated in complete DMEM media at 37°C and 5% CO₂. During passaging, media was aspirated and cells were washed once with 2 mL Gibco™ pH 7.5 phosphate buffered saline (PBS) (Thermo Fisher Scientific, USA, 10010023). Once PBS was aspirated, cells were detached from the tissue culture plastic with 3 mL of HyClone™ Trypsin-EDTA (1X, 0.05%) (Cytiva, USA, SH30236.01) and incubated at 37°C for 5 min. The trypsin was then neutralized with complete DMEM media in a 1:3 ratio. Number of cells was determined using the Beckman Coulter CytoFLEX Flow Cytometer. In each passage, cells were reseeded at a density of 1.7 - 2 x

10^4 cells/cm². Cell seeding density was optimized depending on doubling time, with the ideal doubling time to be 22-24 hours. WT A549s were passaged every 2-3 days with a maximum of 20 passages before bleaching, and were frozen down within the first 2 passages.

To freeze down A549 cells, 1×10^6 cells were pelleted and then added in freezing media, composed of FBS and 10% v/v Anhydrous DMSO (Thermo Fisher Scientific, USA, D12345). Cells were frozen at -80°C for 24 hours, then moved to a liquid nitrogen tank (-160°C).

PolyHEMA-Covered Plates

To make polyHEMA stock, 1.3 g of poly(2-hydroxyethyl methacrylate) (poly-HEMA) (SigmaAldrich, USA, P3932) was added to 33 mL of 99% ethanol at 37°C and left to be stirred overnight until fully dissolved, which was then stored at room temperature. To make polyHEMA-covered plates, polyHEMA stock was added to tissue culture plates (50 μ L for each well in a 96-well plate, 200 μ L for each well in a 24-well plate, 800 μ L to each well in a 6-well plate) and dried overnight in room temperature. Plates were UV-sterilized for 30 minutes and stored in 4°C for up to 1 month before use.

Spheroid Formation

Spheroid formation protocol was adapted from Han *et al*, 2020 [7]. Cells were seeded in polyHEMA-coated tissue culture plates Cells were seeded at two different densities, 50,000 (or 50k) cells/cm² and 75,000 (or 75k) cells/cm² in complete DMEM media with three different methylcellulose (SigmaAldrich, USA, M0512) concentrations, 0.75%, 1%, and 1.25%. Total cell-media mixture volume was 500 μ L/cm² for each well. Methylcellulose prevents excessive aggregation of cells in spheroid culture and maintains even spheroid size. Spheroids were formed

after 3-5 days of incubation, and pictures were taken using BioTek Cytation and analyzed on ImageJ for spheroid size and count. Cell viability was determined using CellTiter-Glo® Luminescent Cell Viability Assay (Promega, USA, G7570) following manufacturer's protocol and ATP luminescence was measured on Biotek Microplate Reader. Optimal condition was determined based on high spheroid size and count as well as cell viability.

For subsequent experiments, cells were seeded at 75k cells/cm² in complete DMEM with 0.75% methylcellulose on 6-well plates. Spheroid doubling time was determined after 5 days. Cells were dissociated into single cells Invitrogen™ Accutase (Thermo Fisher Scientific, USA, 00-4555-56) and counted. Accutase-dissociated spheroids were used for RT-qPCR and γ H2AX staining.

Spheroid Accutase Dissociation

5 days after cell seeding, spheroids were collected into different tubes and plates were washed twice with PBS. Then, additional PBS was added so that original cell culture media was diluted 1:3. Cells were collected at 800g for 15 minutes at room temperature. Cell pellet was resuspended in 1 mL Accutase, then incubated 37°C and 5% CO₂ for 10 minutes. 3 mL of complete DMEM was then added to neutralize the Accutase.

γ H2AX Staining

Spheroids were formed in 6-well plates and 2D A549s were seeded at 1-2 million cells/well. Half the cell population was induced with 4 μ M of etoposide and all four populations of cells (2D cells with and without drug treatment and 3D cells with or without drug treatment) were incubated at 37°C and 5% CO₂ for at least 24 hours. Spheroids were dissociated, and alongside 2D cells, were then fixed at 37°C for 10 minutes. Cells were washed with 500 μ L of

blocking buffer and then permeabilized on ice for 30 minutes. After permeabilizing, cells were then washed with blocking buffer before a 1:1000 antibody dilution was added. Samples were shaken at room temperature for 60 minutes then washed with blocking buffer. Pellets were resuspended in 150 μ L of blocking buffer for flow cytometry.

RNA Extraction and SOX2 RT-qPCR

Spheroids were dissociated, and their RNA along with RNA of 2D cells were extracted using QIAGEN RNeasy Mini Kit (74014) following the manufacturer's protocol. RNA amount was quantified using NanoDrop One Microvolume. 1 μ g of RNA was transcribed into cDNA following Hess Lab protocol using Oligo(dT)₁₈ Primer (ThermoFisher, USA, SO132), 5X AMV RT Buffer (Promega, USA, M5101), AMV Reverse Transcriptase (Promega, USA, M5101), 10 mM dNTP Mix (ThermoFisher, USA, R0192), and RNasin® Ribonuclease Inhibitor (Promega, USA, N2111). After cDNA synthesis, 50 μ L of H₂O was added to the cDNA mix. qPCR was performed on 1:4 diluted cDNA solution following TaqMan® Gene Expression Assays protocol in TaqMan® Array Standard 96-well Plates (20 μ L Reaction) (ThermoFisher, USA, 4391524), with thermal cycling conditions following Applied Biosystems 7500 Fast Real-Time PCR System, standard run, 96-well) [8]. Primers used were SOX2-OT (FAM reporter) (ThermoFisher, USA, Assay ID Hs00415716_m1) and GAPDH (VIC reporter) (ThermoFisher, USA, 4310884E).

Results

Spheroid Formation

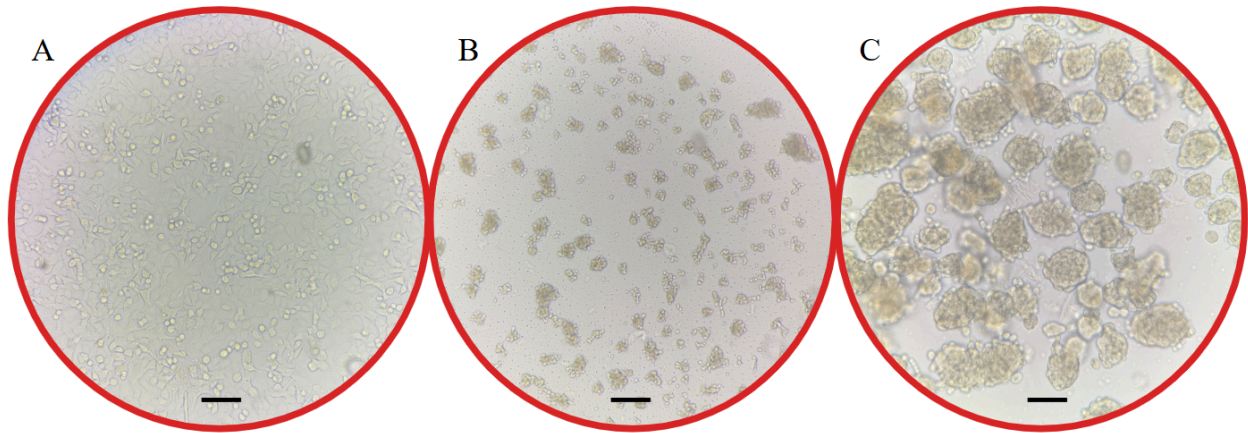


Figure 1. A) WT A549 cells at 147% confluence under 20x, brightfield microscopy, B) WT A549 seeded at 75k cells/cm² spheroids under 10x, brightfield microscopy. C) WT A549 spheroids at day 5 under 20x, brightfield microscopy. Scale bars = 100 μ m.

A549 cells were seeded at varying cell densities and methylcellulose concentrations to determine the optimal condition for subsequent experiments. The ideal condition yielded relatively high spheroid size and spheroid count, as well as high cell viability.

While the average spheroid diameters among conditions were similar, more noticeable differences were seen in the average spheroid count (Figure 2B). Significant differences in spheroid diameter were observed between groups seeded at different cell densities and methylcellulose concentrations. A cell density of 75k cells/cm² and methylcellulose concentration of 0.75% yielded the highest mean spheroid diameter of 89.4 μ m (Figure 2A) and mean spheroid count of 124 (Figure 2B).

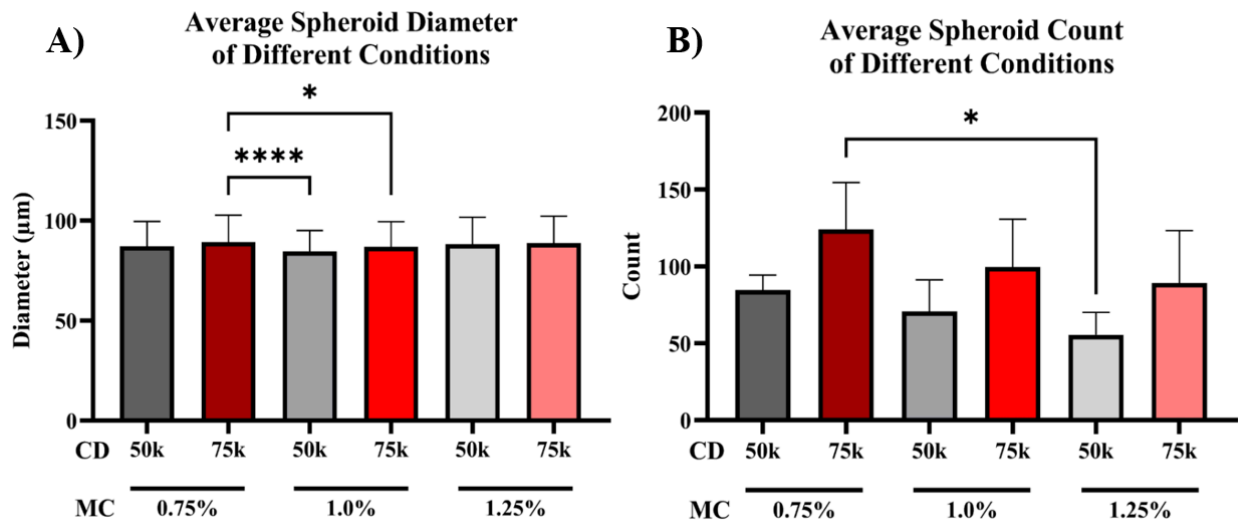


Figure 2. A) Average spheroid size for each condition (n=4). B) Average spheroid count for each condition (n=4). All data are presented as mean \pm std. CD: cell density (cells/cm²), MC: methylcellulose concentration, *: p < 0.05, ****: p < 0.0001. Comparison with condition 75k cells/cm², 0.75% methylcellulose. Non statistically significant, not shown.

To determine cell viability, CellTiter-Glo® Luminescent Cell Viability Assay, which shows ATP luminescence positively correlating with live cell count, was conducted. Average luminescence results from the assay of each condition were normalized by subtraction from the average luminescence of just media. The 75k cells/cm² seeding density and 0.75% methylcellulose concentration condition had the highest luminescence, indicating the condition had the highest cell viability (Figure 3).

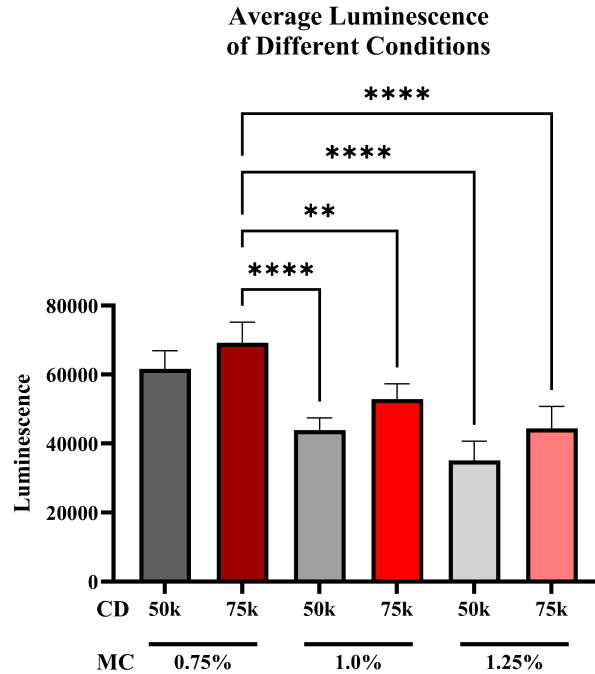


Figure 3. Average luminescence from CellTiter-Glo® Luminescent Cell Viability Assay for each condition (n=4). All data are presented as mean \pm std. CD: cell density (cells/cm³), MC: methylcellulose concentration, **: p < 0.01, ****: p < 0.0001. Comparison with condition 75k cells/cm², 0.75% methylcellulose if not shown is not statistically significant.

SOX2 RT-qPCR

To confirm that the formed spheroids adequately mimic *in vivo* environments, genetic expression of spheroids was analyzed using reverse transcriptase quantitative polymerase chain reaction (RT-qPCR). SOX2 was selected as the gene of interest due to its established association with stemness and tumor progression. Specifically, SOX2 expression has been found to be upregulated in A549 spheroid cells [3]. Further, GAPDH was selected as a housekeeping gene due to its relatively stable expression levels across cell types and culture conditions [9]. RT-qPCR results indicated that SOX2 was not amplified for either cells cultured as single cells (hereafter referenced as 2D) or cells cultured as spheroids (hereafter referenced as 3D) (data not

included). Further, average GAPDH Ct values were not significantly different between 2D and 3D conditions (Figure 4).

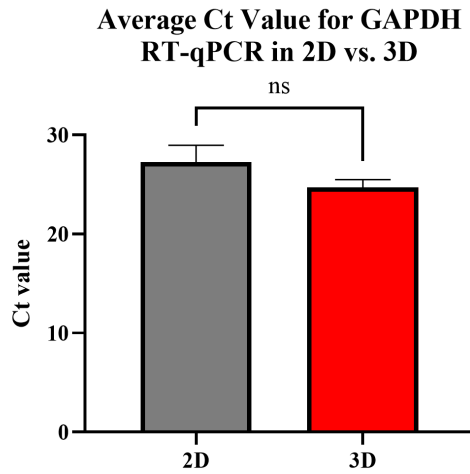


Figure 4. Ct values for GAPDH RT-qPCR for 2D vs. 3D (n=11). Data are presented as mean \pm std. ns: not significant.

γ H2AX Staining

To compare DNA damage caused by Etoposide, 2D WT A549s, Etoposide-treated 2D WT A549s, 3D WT A549s, and Etoposide-treated 3D WT A549s populations were stained to detect DSBs. The fluorescence data gathered from the flow cytometer was then transferred to FlowJo and analyzed. In order to filter out baseline DNA damage ($-\gamma$ H2AX), a threshold value of 10^4 was set and applied to all samples. Any gated values above 10^4 were considered $+\gamma$ H2AX. The luminescence data below in Figure 5 shows similar $-/+ \gamma$ H2AX peak intensities (APC) for 2D+drug (tube 7) and 3D+drug (tube 1). However, when all $+\gamma$ H2AX values for all samples were normalized and compared, there was no statistical significance between the populations (results not shown).

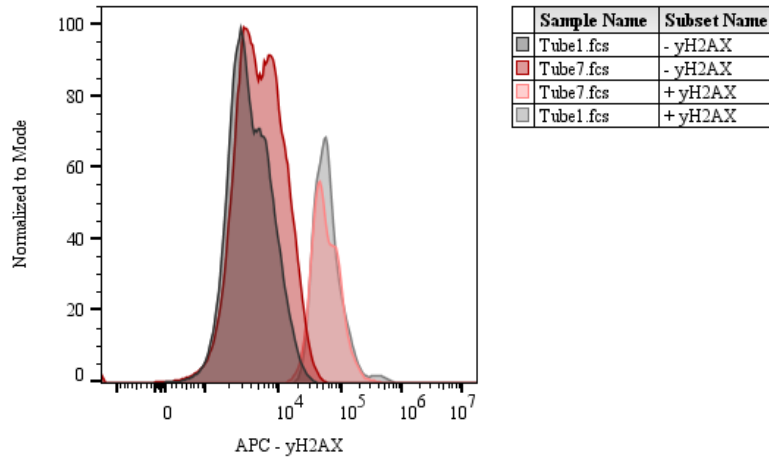


Figure 5. Single-tube APC value comparison for 2D+drug (tube 7) and 3D+drug (tube 1).

The +/-γH2AX APC values gathered from the data gates created on FlowJo were then averaged and a ratio of +γH2AX/-γH2AX was graphed. This normalized data below in Figure 6 showed a statistically significant difference, $p < 0.05$, between the 2D+drug and 3D+drug populations.

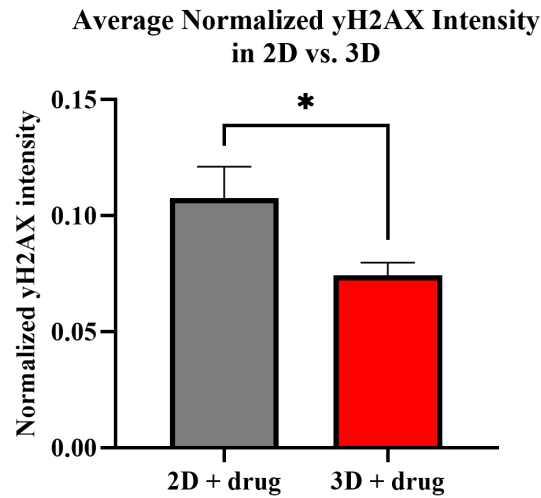


Figure 6. Average normalized APC intensity value comparison for 2D drug (n=3) and 3D drug (n=3). Data are presented as mean \pm std. * $p < 0.05$.

Discussion

To confirm preliminary findings regarding A549 spheroid viability, optimal seeding density, and optimal methylcellulose concentration, spheroids were formed using polyHEMA-coated plates at the six different conditions outlined previously. Average spheroid size was then analyzed for each condition via ImageJ and spheroid cell viability via a CellTiter-Glo® Luminescent Cell Viability Assay. Results indicated that the 75k cells/cm² seeding density and 0.75% methylcellulose condition presented a seeding density and methylcellulose concentration that allowed for optimal spheroid formation. Specifically, spheroids formed under these conditions were significantly greater in spheroid diameter than spheroids formed under 1.0% methylcellulose and either 50k cells/cm² or 75k cells/cm² seeding density (Figure 2A). Spheroids formed at 75k cells/cm² and 0.75% methylcellulose were also significantly greater in number compared to spheroids formed at 50k cells/cm² and 1.25% methylcellulose (Figure 2B). Furthermore, spheroids formed at 75k cells/cm² and 0.75% methylcellulose had significantly greater cell viability than in all but the 50k cells/cm² and 0.75% methylcellulose condition (Figure 3). Based on these findings, the 75k cells/cm² and 0.75% methylcellulose condition was selected to be used moving forward for spheroid formation as it would ensure an adequate spheroid size, count, and viability while conserving methylcellulose and reaching a cell density scale required for subsequent experiments.

Using the optimized parameters above, the spheroid dissociation protocol from Han *et al*, 2020 was optimized to ensure maximum retention of WT A549 cells for subsequent experiments [7]. Specific steps added in conjunction with the rest of this protocol to maximize cell retention included pipetting gently up and down the sides of each well to mechanically break up spheroids and washing each well with PBS. Further, spheroids were dissociated in 1mL of Accutase for 10

minutes, after which a 3:1 ratio of DMEM to Accutase was added. These last two steps were performed to ensure that spheroids were dissociated quickly and completely, while ensuring that the Accutase was neutralized after sufficient time for dissociation. From this dissociation protocol, an estimated spheroid doubling time of 1.5 times per 5 days was determined using a CytoFLEX Flow Cytometer. Of note, it was assumed that all live cells counted from the CytoFLEX were dissociated from spheroids, and that any remaining single cells were counted as dead cells. Though it is possible that this assumption could have affected the true doubling time of spheroids, a sufficient number of cells were obtained for further spheroid culture and experiment from following this protocol.

RT-qPCR results indicated that SOX2 was not amplified on either 2D or 3D conditions. Given that it is well established that Ct values below 29 indicate a high amount of target sequence, results indicated that GAPDH was sufficiently amplified on both 2D and 3D conditions (Figure 4) [9]. Similar GAPDH Ct values in this study indicated that genetic information was efficiently obtained in both conditions for the RT-qPCR assay, and thus the lack of SOX2 amplification was likely not due to errors in RNA extraction or cDNA synthesis.

γ H2AX staining results indicated the Etoposide concentration used was not sufficient to induce DNA damage in 2D and 3D cells (Figure 5). Although these results failed to confirm the hypothesis that there was an increase in DSBs in etoposide versus control and 3D versus 2D, the similar $-/\gamma$ H2AX peak intensities (APC) for 2D+drug versus 3D+drug (Figure 6) validated the success of the staining protocol despite the drug concentration being insufficient to induce a statistically significant difference in DNA damage. Of note, A549 cells in passage 19 and 20 were used for staining 2D and 3D cells respectively, which could have impacted results observed

due to characterized alterations and further deviation from original tumors seen with cells at high passage numbers [10].

Future work will consist of troubleshooting the SOX2 RT-qPCR assay to ensure it is working efficiently under this study's conditions and cell line. This could involve purchasing the assay from a separate site that ensures primer and probe sets have been tested on selected cell lines. Antibody staining for SOX2 could also be performed as an alternative option to ensure SOX2 is being expressed in these A549 spheroids. Further investigation is also needed for utilizing a higher etoposide concentration for γ H2AX staining to ensure efficient induction of DNA damage. Specifically, concentrations of 10 μ M and 100 μ M would be utilized, as these have been previously validated to effectively induce γ H2AX-detectable DNA damage in A549 cells [11]. If such etoposide concentrations do not work, a different chemotherapeutic agent like gemcitabine or a lower passage number may be utilized for the γ H2AX staining protocol to seek clearer results. Finally, future investigations with γ H2AX staining will be performed with an A549 CRISPRi cell line to prepare for genome-wide screening [12]. Ultimately, by employing γ H2AX staining in a 3D environment, more accurate predictions could be made about genes contributing to tumor progression and metastasis.

Acknowledgments

A special thank you to Ms. Carley Schwartz, Dr. Gaelen Hess, and Dr. Paul Campagnola for their contributions to our research. All funding was provided by the Hess Laboratory, located in the Wisconsin Institute of Medical Research (WIMR) at the University of Wisconsin-Madison.

References

- [1] “Lung Cancer Survival Rates | 5-Year Survival Rates for Lung Cancer.” Accessed: Feb. 28, 2025. [Online]. Available: <https://www.cancer.org/cancer/types/lung-cancer/detection-diagnosis-staging/survival-rates.html>
- [2] J. M. Rozenberg *et al.*, “Biomedical Applications of Non-Small Cell Lung Cancer Spheroids,” *Front. Oncol.*, vol. 11, p. 791069, Dec. 2021, doi: 10.3389/fonc.2021.791069.
- [3] C. Choe, H. Kim, S. Min, S. Park, J. Seo, and S. Roh, “SOX2, a stemness gene, induces progression of NSCLC A549 cells toward anchorage-independent growth and chemoresistance to vinblastine,” *OncoTargets Ther.*, vol. 11, pp. 6197–6207, Sep. 2018, doi: 10.2147/OTT.S175810.
- [4] W. M. Bonner *et al.*, “ γ H2AX and cancer,” *Nat. Rev. Cancer*, vol. 8, no. 12, pp. 957–967, Dec. 2008, doi: 10.1038/nrc2523.
- [5] S. Riffle, R. N. Pandey, M. Albert, and R. S. Hegde, “Linking hypoxia, DNA damage and proliferation in multicellular tumor spheroids,” *BMC Cancer*, vol. 17, p. 338, May 2017, doi: 10.1186/s12885-017-3319-0.
- [6] C. Luo, P. J. Johnston, S. H. MacPhail, J. P. Banáth, A. Oloumi, and P. L. Olive, “Cell fusion studies to examine the mechanism for etoposide resistance in Chinese hamster V79 spheroids,” *Exp. Cell Res.*, vol. 243, no. 2, pp. 282–289, Sep. 1998, doi: 10.1006/excr.1998.4170.
- [7] K. Han *et al.*, “CRISPR screens in cancer spheroids identify 3D growth-specific vulnerabilities,” *Nature*, vol. 580, no. 7801, pp. 136–141, Apr. 2020, doi: 10.1038/s41586-020-2099-x.
- [8] “TaqMan Array Plates - US.” Accessed: May 05, 2025. [Online]. Available: <https://www.thermofisher.com/us/en/home/life-science/pcr/real-time-pcr/real-time-pcr-assays/taqman-gene-expression/96-well-taqman-gene-expression-assays.html>
- [9] “GAPDH as a housekeeping gene: analysis of GAPDH mRNA expression in a panel of 72 human tissues | Physiological Genomics | American Physiological Society.” Accessed: May 05, 2025. [Online]. Available: <https://journals.physiology.org/doi/full/10.1152/physiolgenomics.00025.2005>
- [10] “Passage number of cancer cell lines: Importance, intricacies, and way-forward - Prasad - 2023 - Biotechnology and Bioengineering - Wiley Online Library.” Accessed: May 05, 2025. [Online]. Available: <https://analyticalsciencejournals.onlinelibrary.wiley.com/doi/abs/10.1002/bit.28496>
- [11] F. K. Noubissi, A. A. McBride, H. G. Leppert, L. J. Millet, X. Wang, and S. M. Davern, “Detection and quantification of γ -H2AX using a dissociation enhanced lanthanide fluorescence immunoassay,” *Sci. Rep.*, vol. 11, no. 1, p. 8945, Apr. 2021, doi: 10.1038/s41598-021-88296-3.
- [12] T. Fu *et al.*, “p21 promotes gemcitabine tolerance in A549 cells by inhibiting DNA damage and altering the cell cycle,” *Oncol. Lett.*, vol. 26, no. 5, p. 471, Sep. 2023, doi: 10.3892/ol.2023.14059.
- [13] “Imager, Reader Control, & Analysis Software | Agilent.” Accessed: [Online]. Available: <https://www.agilent.com/en/product/microplate-instrumentation/microplate-instrumentation-control-analysis-software/imager-reader-control-analysis-software>

Detailed Methods

Spheroid Formation Protocol

Cell Thawing

1. Prepare 10mL of the appropriate media in a 15ml tube.
2. Using a P1000 slowly mix the media in the cryovial to thaw the pellet of cells.
3. Set the pipette to something below 1000uL, maybe 900uL to avoid clogging the filter while resuspending
4. Once thawed, add to the remaining media and spin down at 1400 rpm for 5 min.
5. Resuspend the cells in fresh media and add to a flask to incubate for typically 3 days before passaging further.

Passaging Adherent (eg. WT A549) Cell Line

1. This assumes passaging in a T75 flask, if using a T150, double the volumes)
2. Remove all media with Pasteur pipette and vacuum
3. Add 2ml of PBS to remove remaining media and remove with vacuum
4. Add 2mL trypsin and wait five min at 37C
5. Add 6mL media to neutralize trypsin while mixing to remove cells from the bottom of the flask (at least 3x the volume of trypsin used)
6. Wash the flask with the cells in media, remove and add to a 15ml conical tube
7. Return enough media to the flask so cells will be confluent when next checked
8. Add additional media (eg. 10 ml)
9. The media just needs to be enough to cover the cells, provide nutrients and won't dry out. Adding more media doesn't change the confluence of the cells.

Counting Cells

1. Add 100uL of cells in media to a 1.5 mL eppendorf tube
2. Measure 10uL on the cytometer and gate for the live population of cells
3. Use the "live" number of events to calculate the confluence and total number of cells you have.
4. Since the cells have been diluted 1:10, the "live" events gives you the number of cells per uL, multiply this by 1000 to get the number of cells per mL.

CytoFLEX Cytometer

Once the CytoFLEX and its corresponding computer have been turned on, the following steps must be followed:

- 1) In the top left click Cytometer > Daily clean
 - a) Load Blue tube and press run (3 min)
 - b) Load DI water and press run (1 min)
- 2) To create a new experiment click file> new> from template>Browse
 - a) Click into BME Team cell counting file on computer drive Documents>CytoFLEX>BME Team
- 3) Name file after current date and passage number
- 4) Vortex sample
- 5) Uncap sample in 1 mL Eppendorf container and place in CytoFLEX
- 6) Select fast in settings on left side
- 7) Click record
- 8) Once a dot plot of the events has been created right click on the plot>properties
 - a) Auto set x-axis
 - b) Auto set y-axis
- 9) Right click on graph and select polygon tool
- 10) Gate the cells in a oblong hexagon fashion trying to bin all relevant events (general note: live cells are larger than dead cells)
- 11) Save file
- 12) Rerun daily clean (step 1)

Methylcellulose Stock Solution

Materials:

- 2 g Methylcellulose
- 100 mL deionized water (diH₂O)
- 1x beaker that can hold 100 mL
- 1x heated stirrer plate
- 1x magnetic stir bar
- 1x 500 mL glass bottle with screw cap

Steps:

1. Add stir bar into glass bottle
2. Add methylcellulose and diH₂O into glass bottle
3. Mix solution at 80°C until homogeneous
4. Autoclave solution using Liquid 2 (L2) setting with screw cap loose
5. Move solution into walk-in fridge and mix overnight at 4°C
6. Keep solution at 4°C until needed

Poly-HEMA Stock Solution Preparation

Materials

- 1.3 g poly-HEMA
- 33 mL 99% ethanol

- 1x beaker that can hold 33 mL
- 1x heated stirrer plate
- 1x magnetic stir bar
- 1x 0.2 μm filter (Nalgene) or something similar
- 1x 250 mL glass bottle with screw cap

Steps:

1. Add poly-HEMA to 99% ethanol
2. Mix solution overnight at 37°C
3. Filter solution into glass bottle
4. Keep solution at 20°C (room temperature) until needed

Cell Seeding

Materials:

- Methylcellulose stock solution
- Poly-HEMA stock solution
- Cells collected after passaging
- Serum-free DMEM
- 15 mL conical tubes
- Eppendorf tubes
- 96-well plates

Procedure:

- 1) Collect cells from that day's passage and obtain cell concentration using the CytoFLEX to obtain cell concentration (cells/mL).
- 2) Obtain a 1.5 mL cell solution of 500,000 cells/mL (750,000 cells in total)
 - a) Volume needed to obtain 750,000 cells: $750,000 \text{ (cells)} / \text{cell concentration from CytoFLEX (Step 1, cells/mL)}$
 - b) Obtain the volume calculated from the T75 flask and transfer into a 15 mL conical tube.
 - c) Spin down the the tube at 200g for 5 min, then remove supernatant to collect cells using vacuum filter
 - d) Resuspend cells in 1.5 mL of serum-free DMEM.
- 3) Prepare master mix of cell solutions for spheroid formation for 4 different densities. The master mix can be used for 6 wells in a 96-well plate and can be scaled to adjust for different numbers of wells if needed.
 - a) Density 1: 25,000 cells/cm²
 - i) Into a new 1.5 mL eppendorf tube, add:
 - (1) 100 μL of 500,000 cells/mL cell solution (Step 2)
 - (2) 525 μL of serum-free DMEM
 - (3) 375 μL of 2% methylcellulose stock
 - ii) Mix gently by pipetting up and down, make sure the solution is well-mixed
 - b) Density 2: 50,000 cells/cm²
 - i) Into a new 1.5 mL eppendorf tube, add:

- (1) 200 μL of 500,000 cells/mL cell solution (Step 2)
 - (2) 425 μL of serum-free DMEM
 - (3) 375 μL of 2% methylcellulose stock
 - ii) Mix gently by pipetting up and down, make sure the solution is well-mixed
 - c) Density 3: 75,000 cells/cm²
 - i) Into a new 1.5 mL eppendorf tube, add:
 - (1) 300 μL of 500,000 cells/mL cell solution (Step 2)
 - (2) 325 μL of serum-free DMEM
 - (3) 375 μL of 2% methylcellulose stock
 - ii) Mix gently by pipetting up and down, make sure the solution is well-mixed
 - d) Density 4: 150,000 cells/cm²
 - i) Into a new 1.5 mL eppendorf tube, add:
 - (1) 600 μL of 500,000 cells/mL cell solution (Step 2)
 - (2) 25 μL of serum-free DMEM
 - (3) 375 μL of 2% methylcellulose stock
 - ii) Mix gently by pipetting up and down, make sure the solution is well-mixed
- 4) Aliquot the prepared master mix into each well: add 160 μL of cell mixture into each well.
- 5) Grow spheroids in 37C, 5% CO₂ incubator. Spheroids are split every 3-4 days.

References

[1] K. Han *et al.*, “CRISPR screens in cancer spheroids identify 3D growth-specific vulnerabilities,” *Nature*, vol. 580, no. 7801, pp. 136–141, Apr. 2020, doi: 10.1038/s41586-020-2099-x.

BioTek Cytation Imaging

Started BioTek Gen5 3.14 application, placed 96 well plate in BioTek with A1 well on the bottom right of the holder, created new experiment and new protocol— made sure that the experiments and photos are saved in the same location and transferred together— set protocol and imaging settings, and clicked run [13]. The protocol procedure was set as follows: temperature at 37°C; middle imaging; deselected auto exposure; selected plate type; and set imaging settings to 4x magnification, Brightfield, FVOW, set z-stack, and distance. Once the experiment trial was run and images and experiments were saved together in the same file, a scale bar was added.

ImageJ Analysis of BioTek Cytation Images

Out of 120 images taken, 5 z-stack slices per 6 wells for 4 conditions, take the most in-focus images ("Spheroid Images" entry) complete the following steps:

1. Open ImageJ application and upload desired image

2. Set scale bar using the scale bar on the image
3. Go to image → type → 16 bit
4. Threshold image
5. Go to process → binary → watershed
6. Trace around well circle to isolate spheroids as the only thing being analyzed
7. Go to analyze particles → set interval to 4147-infinity, select show outlines from drop down, check “display results” and “exclude on edges” boxes, clear results boxes, and click OK
 - a. 4147 μm^2 taken from measurement of spheroid with ≥ 20 cells which is deemed the benchmark to accurately model the tumor micro environment [1].
 - b. See entry "Spheroid Size Threshold_ImageJ Analysis_11/25/2024" in Emily Rhine folder
8. Save measurements to an excel sheet
9. Record data results to separate excel file by combine measurement averages
10. In excel compiled sheet and run 4 way ANOVA test in excel
11. Generate a suitable plot for this comparison (see example below)

CellTiter-Glo® luminescent cell viability assay

Follow Promega CellTiter-Glo® luminescent cell viability assay instruction manual:

1. Add 1X volume of Cell Titer Glo mix (stored at -30C) to well volume.
2. Mix plates for 2 minutes at RT on an orbital shaker. Cover plate to reduce light exposure.
3. Incubate plates for 10 minutes at room temperature. Cover plate to reduce light exposure.
4. This plate can now be analyzed on a chemiluminescence machine (peak emission wavelength of cell titer glo is 560nm).
5. Analyze data: Normalized triplicate cell viability numbers for each drug treatment to vehicle only numbers (column 2). Plot on a log-linear plot and determine approximate IC50 values.

Spheroid Dissociation

1. Prepare 6, 50 mL, conical tubes
2. Pipette contents of wells up and down using P1000 set to 800uL to break up spheroids and move them away from well walls (do this alongside step 3 for each well)
 - a. Note: Try to make as few bubbles as possible in this step
3. Transfer the media (with cancer spheroids inside) from the wells into the conical 50 mL tubes using a P1000 pipet tip set to 1000 μL .
4. Wash well with 1000 uL PBS in each well *twice* to collect any remain cells and add to the corresponding conical tube
 - a. Note: Try to make as few bubbles as possible in this step/
5. To each 50 mL conical tube, add 14.4 mL of PBS
6. Pellet the spheroids via centrifugation at 800g, 15 min, then remove the supernatant using a P1000 pipet

7. Remove all liquid from conical tube 50 mL at a time using P1000 pipette making sure not to remove the cell pellet from the conical tube
 - a. Note: Try to make as few bubbles as possible in this step
 - b. Note: You don't have to change pipette tips as you remove the media unless you think the tip may be contaminated
 - c. Note: do this step concurrently with step 6 (the step after)
8. Resuspend the pellets in 1000 μ L Accutase using a P1000 pipette.
9. Thoroughly mix the solution by pipetting up and down 5 times and going against the wall of the tube (or until you can visibly see the cell pellet has been broken up a little bit) using a P1000 pipette
 - a. Note: Complete this step concurrently with the last step so that the cells do not dry out between the steps. Once one tube is aspirated add accutase right away.
10. Incubate the tubes in 37°C for 10 minutes.
11. Pipette the solution in each tube up and down 5 times using a P1000 pipette set to 1000 μ L, add 3000 μ L of DMEM so that the final volume reaches 4 mL

γ H2AX Stain

Materials:

- Fix buffer (Fisher Scientific, BDB557870)
- Permeabilization buffer III (Fisher Scientific, BDB558050)
- Blocking buffer (10% FBS in PBS)
- Primary conjugated γ H2AX antibody
- Etoposide

Drug	Solvent	Stock []	Working []
Etoposide	DMSO	16.991mM	4 uM

Day 1 Procedure:

1. Thaw Etoposide
2. Seed 6 well plate well plate with 2 million 2D cells per well
3. Spin down volume of 2D cells equal to 6-12 million cells total in 15 mL conical tube
4. Resuspend cells in 6 mL complete DMEM
5. Transfer 1 mL of cells to each 2D well
6. Add 5 mL of complete DMEM to control wells
7. Create Master mix of drug for 2D and 3D
8. master mix (2D):
9. 15 mL complete DMEM + 4.2 uL etoposide
10. master mix (3D):

11. 3.6 mL complete DMEM + 4.2 uL etoposide
12. Add 5 mL of master mix (2D) to 2D drug wells
13. Add 1.2 mL of master mix to 3D drug wells
14. Let incubate at 37C for at least 24 hours

	Stock Concentration (mM)	Desired Treatment Concentration (uM)	Individual Well Volume (uL)	Volume of Stock Drug Needed (uL)	# of Wells	Volume of Drug for Master Mix (uL)	Volume of Media Already in Well (mL)	Volume of Media for Master Mix (mL)
2D	16.991	4	6000	1.41	3	4.24	14.4	3.6
3D	16.991	4	6000	1.41	3	4.24	14.4	3.6

Day 2 Procedure:

1. Preheat fix buffer I to 37 degrees Celsius. Cool perm buffer on ice. (need to locate these and know where to heat up fix buffer)
2. In TC, prepare cells as though you were planning to dissociate them
3. 2D:
 - a. Aspirate off old media making sure not to touch the bottom with the pasture pipette
 - b. Wash cells with 1 mL of PBS/ well
 - c. Aspirate PBS
 - d. Add 1 mL of Trypsin/ well
 - e. Incubate for 5 minutes at 37C
 - f. Add 3 mL of complete DMEM to neutralize trypsin
 - g. Transfer each well to separate 15 mL conical tube
4. For 3D see Accutase dissociation protocol
5. Centrifuge 15 mL conical tubes at 500 x g for 5 min.
6. Aspirate excess liquid, being careful to avoid disturbing the pellet.

****all following steps will be done at your bench in non-sterile conditions not in the BSC****

7. Resuspend dissociated cells in fix buffer at a ratio of 10 uL fix buffer per million cells. Fix at 37 degrees Celsius for 10-15 min.
 - a. 30 uL/ sample
8. Wash cells with a 500 uL of blocking buffer.

- a. *Note: For all washes, centrifuge at 500 g for 5 min, then aspirate excess liquid, being careful to avoid disturbing the pellet.***
9. Permeabilize by adding 20 uL/million cells of permeabilization buffer to the pellet. Pipet up and down a few times to ensure the pellet is evenly mixed. Permeabilize on ice for 30 minutes.
10. 60 uL/ sample
11. Meanwhile, prepare antibody by making a 1:1000 dilution of antibody in blocking buffer (depending on the brand/type of primary antibody used, this ratio may differ).
 - a. Math: Antibody solution: 12 samples total, go for overage so make for 15. We do 100 uL per every 10^7 cells present which is 10 million cells. Since we have 2-3 million total, we will do 40uL for each sample.
 - b. 15 samples x 40 uL = 600 uL of blocking buffer
 - c. To do a 1000x dilution: $600/1000 = 0.6$ uL of antibody
12. After permeabilizing, wash the samples once with blocking buffer.
13. Stain with the antibody dilution at an amount of 100 uL mixture per 10 million cells. Place on a Nutator (shaker plate) at room temperature for 1 hr.
 - a. $(100\text{uL}/10 \text{ million}) \times (\text{xuL}/3\text{million}) = 30\text{uL}/ \text{sample tube}$
14. Wash with blocking buffer (500 uL of blocking buffer, spin down, and aspirate)
15. Resuspend pellet in 150 uL of blocking buffer for flow cytometry.

Fluorescence Reading using Flow Cytometry

1. Set up fluorescence reading
2. Run QC
3. Settings:
4. Open new experiment and name appropriately
5. Alter stop to **150,000** events
6. Change to fast read
7. Change stop volume to 140uL
8. Run daily clean
9. Load first sample (3D drug - most cell death expected)
10. If after 2 minutes less than 150,000 events have been read, change events to match expected events for the rest of the samples **making sure not to set it less than 50,000 events.**
11. Run the rest of the samples.
12. Save data to Hess Lab drive
13. Run daily clean

RT-qPCR

Step 1: RNA extraction

1. Clean bench and spray down with RNaseZap
2. Make RLT buffer (350µL per sample) by adding β-mercaptoethanol (10µL for every 1 mL of RLT buffer) into a separate tube inside the fume hood.

3. Note: β -mercaptoethanol should be added freshly before each set of preparations.
4. Prepare 1-5 million cells and spin down at 500g for 5 minutes. Aspirate supernatant and resuspend cells in 350 μ L of prepared RLT buffer.
5. Pipet cell solution onto a QiaShredder column and spin for 2 minutes at 14000 rpm. Discard column, not flow-through!
6. Add 350 μ L of 70% ethanol to flowthrough. Transfer mixture to RNAeasy spin column and spin for 30s at 10000rpm. Discard flow-through.
7. Prepare Dnase solutions by adding 10 μ L of DNase I stock solution to 70 μ L of RDD buffer for each column. Mix well.
8. DNase in the resuspended form is stored in the -30 in the box with reagents for cDNA synthesis
9. Add 80 μ L of DNase I + RDD mix onto the column and let sit at room temperature for 15 min.
10. Add 350 μ L of RW1 buffer to column and spin 30s at 10000rpm. Discard flowthrough.
11. Add 500 μ L of RPE buffer to column and spin 30s at 10000rpm. Discard flowthrough.
12. Add 500 μ L of RPE buffer to column and spin 2 minutes at 14000rpm. Discard flowthrough.
13. Place Column into new 1.5mL tube and elute with 30-100 μ L of RNase-free water. Let sit for >5 minutes and spin down for 1 minute at 10000 rpm. Nanodrop to determine concentration.
14. Purified RNA can be stored at -80 C for 6 months, -20 C for 1 month.

Step 2: cDNA synthesis

1. This protocol is for 20 μ L/1 μ g of RNA protocol. Prepare in PCR tubes:
 - a. RNA (1 μ g) x
 - b. Oligo dT (0.5 μ g/ μ L) 2
 - c. H₂O up to 10 μ L
2. Heat sample at 70°C for 5 minutes. Place heated sample on ice for 5 min.
3. Add 10 μ L of following mix:
 - a. 5X AMV RT Buffer 4
 - b. 10mM dNTPs 1
 - c. RNAsin 0.25
 - d. AMV RT 0.75
 - e. H₂O 4
4. Incubate at 42°C for 2 hours.
5. If going into a PCR, add 50 uL of Milli-Q-H₂O and proceed to set up qPCR reaction

Step 3: qPCR

1. Dilute cDNA solution from Step 3 1:3 by adding 10 uL of cDNA solution into 30 uL of molecular grade H₂O
2. Set up Master Mix 1 (everything but cDNA, 35X):
 - a. 35 uL of 20X Taq Man Gene Expression Assay Target Gene (SOX2)
 - b. 35 uL of 20X Taq Man Gene Expression Assay Housekeeping Gene (GADPH)
 - c. 350 uL of 2x Taq Gene Expression Master Mix

- d. 140 uL of RNase Free water
3. Aliquot 80 uL of Master Mix 1 into each PCR tube. We will have 6 tubes in total pertaining to 6 samples. This is to make Master Mix 2
4. In each tube, add 20 uL of cDNA (diluted 1:4) relative to each sample. We will have a total of 100 uL of Master Mix 2, pertaining to 5 replicates/reactions per sample (5X of 20uL/reaction)
5. Aliquot 20uL of Master Mix 2 into each well
6. Setup the following protocol in the thermal cycler:
 - a. 95°C 10min
 - b. 95°C 15s
 - c. 60°C 60s
 - d. plate read
 - e. Go to b 39 times

Appendix

A. γ H2AX Data from FlowJo

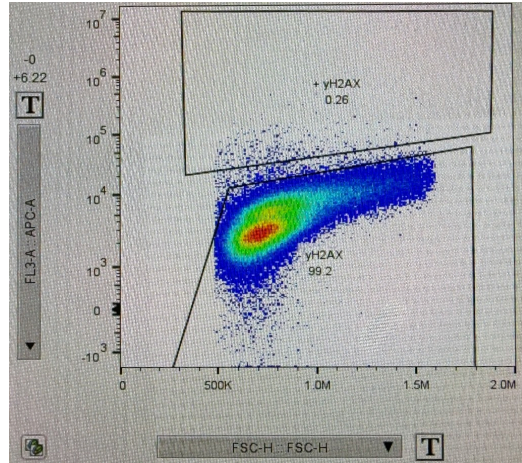


Figure 1. Tube 1 flow cytometry data gaited on FlowJo. Please note the <1% +yH2AX.

Table 1. Raw yH2AX APC Intensity data.

	Alive/+ yH2AX Mean (FL3-A :: APC-A)	Alive/- yH2AX Mean (FL3-A :: APC-A)	yH2AX/-yH2AX
Tube1.fcs	71667	5040	0.0703252543
Tube2.fcs	60891	4906	0.08057019921
Tube3.fcs	62261	4452	0.07150543679
Tube4.fcs	66149	4638	0.07011443862
Tube5.fcs	57622	4830	0.08382215126
Tube6.fcs	60008	4398	0.07329022797
Tube7.fcs	63608	7218	0.1134762923
Tube8.fcs	69854	6418	0.09187734417
Tube9.fcs	61576	7208	0.1170585943
Tube10.fcs	179503	6849	0.03815535116
Tube11.fcs	71783	6876	0.09578869649
Tube12.fcs	74016	7862	0.1062202767

Table 2. Normalized yH2AX APC Intensity data.

	3D drug	2D drug
Mean	0.0741336301	0.1074707436
SD	0.005605379032	0.0136225475

B. Product Design Specifications

CRISPRi Screening in Cancer Spheroids to Investigate Factors in Genome Stability

Updated: January 31, 2025

Function: Although previous CRISPRi screening in 2D monolayers has provided useful knowledge on cancer drivers and therapeutic susceptibilities, it lacks an element of biological relevance to an *in vivo* environment. Therefore, our team was tasked with developing a cell culture method that is compatible with a 3D environment and CRISPRi screening in order to identify sources of DNA mutations in the tumor environment. On a high scale, the team must use the selected cell line, A549, to optimize a spheroid formation protocol, and develop a protocol to stain for γ H2AX: a histone variant that acts as a sensitive marker for DNA damage.

Client requirements:

1. The team must choose an appropriate human cell line for the project, ensuring to maximize cell spheroid viability.
2. The team must develop a spheroid formation protocol for the chosen cell type or line that optimizes seeding density and viscosity.
3. The team must optimize the cancer spheroid characteristics and their formation protocols for CRISPR screens.
4. By May 2025, the team must develop a protocol to stain for γ H2AX: an indicator for the location of maximal damage due to DNA double-strand breaks (DSBs).
5. By May 2025, the team must perform high throughput genome-wide CRISPR screening to check for the effects of different genes on the amount of γ H2AX detected.
6. Future steps post-design project include performing an analysis to determine the phenotypic differences between cells grown in monolayers (2D) versus cells grown as spheroids (3D).

Design requirements:

1. Physical and Operational Characteristics
 - a. *Performance requirements:*
 - i. Spheroid formation protocol should be scalable, repeatable, and optimized regarding:
 1. Seeding densities to ensure only 1 single guide (sgRNA) will enter each cell.
 2. Viscosity of growth medium (by changing the amount of methylcellulose) to ensure homogenous spheroid size.
 3. Optimizing reagent percentages of growth medium to ensure maximum cell growth.
 - ii. γ H2AX staining protocol should be scalable, repeatable and optimized regarding the antibody used and amount of reagents used.
 - b. *Safety:*
 - i. Biosafety Level 2 (BSL-2) is necessary for work with agents associated with human tissue and disease including the human lung cancer cell line which the team will be using [1]. The following universal precautions for BSL-2 are used to reduce the risk of bloodborne pathogens:

1. Laboratory personnel receive specific training in handling pathogenic agents and are supervised by scientists competent in handling infectious agents and associated procedures [1].
 2. Access to the laboratory is restricted when work is being conducted [1].
 3. All procedures in which infectious aerosols or splashes may be created are conducted in BSCs or other physical containment equipment [1].
- ii. Cultures, tissues, specimens of body fluids, or potentially infectious wastes are to be placed in a container with a cover that prevents leakage during handling, processing, and storage [2].
 - iii. Although the team does not anticipate the need to remove our cells or cancer spheroids from Hess Lab, the following BSL-2 guideline must be followed.
 1. Biological agents prepared for transport must be placed in a secondary leak-proof, unbreakable carrier. Carriers must have the biohazard label, including the identity of the agent, affixed to the outer surface of the transport container [2].
 - iv. If genome editing leads to the inactivation of tumor suppressors or activation of oncogenes in human cells, then oncogenesis becomes a possibility in the event of exposure [3]. This is a concern when:
 1. Genome editing components are delivered by a method that can enter human cells (e.g., viral vectors that can infect human cells).
 2. All components required for genome editing are likely to be introduced in a single exposure event, such as a needle stick or mucosal splash.
 3. Genome editing is designed to target human oncogenes/tumor suppressors or can target human oncogenes/tumor suppressors due to homology.
 - v. A risk assessment and containment preparations for research with lentiviral vectors should consider (1) the nature of the vector system, (2) transgene insert, and (3) type of manipulations involved [4]. Either BL-2 or enhanced BL-2 will be appropriate for these such experiments. The major risks to be considered for research with HIV-1 based lentivirus vectors are (1) potential for generation of replication-competent lentivirus (RCL) and (2) potential for oncogenesis [4]. These risks can be mitigated by the nature of the vector system (and its safety features) or exacerbated by the nature of the transgene insert encoded by the vector [4].

c. *Accuracy and Reliability:*

- i. A starting cell seeding density of 50,000 cells/cm² will be used. Cells at these densities will be monitored for growth and death rates on a live-cell microscope for 60 hours to confirm a ~30% cell death rate during the initial growth phase of the spheroids [5].
- ii. Cell seeding densities must be such that only 1 single guide (sgRNA) enters each cell to allow for analyzable results. Therefore, a total of

100-150 million cells will be required for each genome-wide CRISPR screen [6].

iii. A minimum of 3 γ -H2AX staining experiments will be performed to establish confidence in results.

d. *Life in Service:*

i. Before preparing them for CRISPRi screening, cells will need to be in culture for approximately 2 months to allow for optimization of the spheroid formation protocol, including the spheroids' size and cell quantity. This time period may consist of thawing a new vial of cells for continuing cell culture if the passage number exceeds 25 [7].

ii. After infection with the selected genome-wide library, a minimum of 6 days will be required before the cells can be further utilized [5].

iii. After generating spheroids from these cells, a minimum of 21 days will be required for the genome-wide screening [5].

e. *Shelf Life:*

i. When not in culture or experimental use, cells will be stored using cryogenics for long-term use. Frozen cells will be kept in cryovials with medium containing a cryoprotective agent like DMSO. These vials will be immersed in liquid nitrogen at a temperature below -135°C [8].

f. *Operating Environment:*

i. Cells needed to generate spheroids for this project will be kept in standard incubation conditions: 37°C , 5% CO_2 , and 90–95% relative humidity [9].

ii. Cells in culture will be maintained in a cell culture medium supplemented with fetal bovine serum (FBS) to deliver essential growth factors. The medium will contain phenol red, which will indicate the need to replace it if it turns yellow (acidic) or purple (alkaline) [9].

iii. For preparation for CRISPRi screening, spheroid cells will be exposed to a pooled single-guide RNA (sgRNA) library and CRISPR associated protein 9 (Cas9) via lentiviral vectors [10].

g. *Ergonomics:*

i. Protocols should create a more efficient and accurate experimental system that better mimics *in vivo* tumor conditions, while proving the effects of γ H2AX in DNA damage.

ii. The design of cell culture methods should be optimized by selecting appropriate cell lines and developing correct spheroid formations to enhance the accuracy, relevance, and scalability of CRISPRi screening in a 3D tumor model.

iii. The γ H2AX staining protocol should be optimized for improved precision. Additionally, γ H2AX foci counting should be performed via high-throughput screening for efficiency.

h. *Size:*

- i. Spheroids will be generated in 60 mm tissue culture plates [5].
- ii. Cell seeding density for spheroids will be restricted to 50,000 cells/cm² to ensure appropriate CRISPRi screening [5].
- iii. To adhere to a spheroid size suitable for γ H2AX staining, spheroids will be kept to ~500 μ m in diameter [11].

i. *Weight:*

- i. Although solutions have been presented for single-cell analysis, literature does not cover this aspect for 3D models such as spheroids [12]. Tools like a flow-apparatus are needed for the accurate measurements of mass density, size and weight of such spheroids [12].

j. *Materials:*

- i. In order to leave room for the engineering design element, the client did not give specific requirements for materials other than cell type to be used for the project.
- ii. The client has recommended that the team select a human lung cancer cell line that (1) adheres well to spheroids and (2) has a high percentage of viable cells.

k. *Appearance:*

- i. Since the spheroids should have an overall appearance that will accurately replicate the 3D tumor microenvironment, the spheroids must be: translucent, porous, uniformly spherical, and have a stiffness analogous to the lung cancer cell line selected [13].

2. Production Characteristics

a. *Quantity:*

- i. By the end of this semester, the team will design, fabricate, and develop two appropriate protocols for the chosen cell line, one for spheroid formation and one for γ H2AX.
- ii. The final protocols should be capable of accurately mimicking the 3D tumor environment for CRISPR screening, including features for effective spheroid formation and reliable γ H2AX staining to assess DNA changes.
- iii. The protocols can be adapted to other similar human cancer cell lines; however, some of the characteristics may be subjected to change.

b. *Target Product Cost:*

- i. The budget for the project is \$1,000. Some of the materials and services that the team expects to use are listed below; however, many are already available in Hess Lab and will not need to be purchased.
 1. Biological materials such as reagents (for various media), antibodies (for γ H2AX staining) and lentivirus for gene insertion.
 2. Obtaining new cell lines if needed.
 3. Designing of CRISPR library.

- ii. This project is a part of an emerging field, so similar products on the market are not available for non-clinical purchasing purposes. Therefore, the target cost of the final prototype is to remain within the budget proposed.

3. Miscellaneous

a. *Standards and Specifications:*

- i. Cell lines used for this project should comply with the requirements outlined in part 1b - *Safety* above. Specifically, cell lines used and their associated reagents and protocols will need to comply with requirements for BSL-2 labs; meaning, they are or will contain bacteria, viruses, and organisms that pose a moderate safety hazard to healthy adult humans [14].
- ii. The use of human cell lines and other organisms in research for this project is outlined in various standards and regulations including Section 29 CFR 1910.1030 (*OSHA Bloodborne Pathogens*), Section 21 CFR 876.5885 (*Tissue culture media for human ex vivo tissue and cell culture processing applications*), and a comprehensive cell line guideline by Geraughty, *et al* (2014):
 - 1. Bloodborne pathogens are microorganisms in the blood that can infect and cause diseases in humans [15].
 - 2. Human cell lines are designated by the FDA as Class II (Special Controls) [16].
 - 3. If applicable, genetic modification processes on the cell line must be detailedly-noted and kept track of, including but not limited to, sequence details, insertion vectors and modes, and antibiotic resistance markers [17]. Following modification using lentivirus, cell lines may be additionally tested to confirm non-infectivity [17].
- iii. Human cell lines that may be needed for this project are already available in the client's lab, or can be obtained from another lab, meaning they already comply with the requirements listed above. In the case that biological research materials are obtained from another lab or institute, an MTA (material-transfer agreements) must be signed by both the provider and recipient to define their rights regarding the materials [17].
- iv. In the unlikely case that new human cell lines must be acquired, additional ethical and authentication considerations must be taken into account:
 - 1. Ethical considerations are outlined in Part 3b - *Patient-related concerns* below.
 - 2. For cell line authentication, short tandem repeat (STR) is the recommended profiling method, outlined in ASN-0002 (*Authentication of Human Cell Lines: Standardization of STR Profiling*) provided by the American National Standards Institute (ANSI) and the American Type Culture Collection (ATCC). STR is a DNA-based sequence method that allows for detection of

mutations and adventitious agents, as well as confirmation of the origin of the interested cell line [18].

b. *Patient-related concerns:*

- i. As cell lines and other tissue materials needed for this project are either already available in the client's lab or in another lab, their ethical and other patient-related considerations have already been considered.
- ii. In the unlikely case that a new cell line is obtained from patients or donors, the Patient Consent Form and associated Patient Information Sheet will be needed to explain the need for the specimen and the purpose and importance of the research [17]. Patient/ donor and original tissue information such as age, sex, clinical history, site of origin, nature of tissue specimen, stage and grade of cancer or other disease/ pathology, and other important information must also be recorded.
- iii. For maximal biological relevance to patient *in vivo* tissue, the team will ensure our spheroids are designed based on critical 3D parameters: stiffness analogous to lung cancer tissue, replication of the complex cell-cell interactions and cell-extracellular matrix (ECM) interactions, and increased accuracy in CRISPRi screening environment to better detect precancerous cells [19], [20].

c. *Competition:*

- i. 3D cancer spheroid models that can be utilized for CRISPR-Cas9 screens have been implemented in previous research, with and without using an exogenous extracellular matrix (ECM).
 1. Takahashi *et al*, 2020 utilized 4% Matrigel™ as a bio-scaffolding to plate A549 and H1437 cells spheroids, which can better model *in-vivo* interactions between cells and the ECM [21]. The study's CRISPRi library gene list as well as its phenotypes have a strong overlap with those from Han *et al*, 2020 (described below), suggesting that inner cells of ECM-based 3D spheroids may be similar to ECM-free 3D spheroids.
 2. Han *et al*, 2020 devised an ECM-free scalable method to propagate 3D spheroids [5]. H23 cells were seeded on pre-treated ultra-low attachment or polyhema-covered plates in RPMI 1640 medium with 0.75% methylcellulose. Similar methods have been found to be implemented in other studies, such as by Ferrarone, *et al* (2024) and by Stiff, *et al* (2024) [22], [23]. A total of around 210,000 single guide RNAs (sgRNAs) were designed for the CRISPRi library with around 10 sgRNAs per gene [5]. Samples were used to calculate growth and/or tested for drug resistance [5].
 3. 3D stem cell models have also been genomically edited via CRISPR-Cas9 to better model various disease phenotypes. Cells can be assembled with or without supporting scaffolds. Some example scaffolds include laminin-rich Matrigel, collagen matrix, or those that are naturally derived from decellularized tissues [24].

- ii. γ H2AX has been used in various research as a marker to monitor DNA damage and repair, such as in understanding the toxicity of cigarette smoke or in understanding the effects of chemotherapeutic treatments in cancer patients [25], [26]. High throughput systems have also been developed to aid with counting γ H2AX foci [26].

References

- [1] CDC, “Biosafety in Microbiological and Biomedical Laboratories—6th Edition,” p. 37, Jun. 2020.
- [2] NIH, “Laboratory Biosafety Manual BSL-2 and BSL 2/3.” Accessed: Sep. 12, 2024. [Online]. Available: <https://ors.od.nih.gov/sr/dohs/Documents/bsl-2-lab-safety-manual.pdf>
- [3] University of Wisconsin-Madison Office of Biological Safety, “Biosafety Guide to CRISPR.” Accessed: Sep. 12, 2024. [Online]. Available: <https://ehs.wisc.edu/wp-content/uploads/sites/1408/2022/04/EHS-BIO-GUI-034-V02.pdf>
- [4] NIH, “Lenti_Containment_Guidance.pdf,” Major RAC Concerns. Accessed: Sep. 15, 2024. [Online]. Available: https://osp.od.nih.gov/wp-content/uploads/Lenti_Containment_Guidance.pdf
- [5] K. Han *et al.*, “CRISPR screens in cancer spheroids identify 3D growth-specific vulnerabilities,” *Nature*, vol. 580, no. 7801, pp. 136–141, Apr. 2020, doi: 10.1038/s41586-020-2099-x.
- [6] V. Chenouard *et al.*, “Excess of guide RNA reduces knockin efficiency and drastically increases on-target large deletions,” *iScience*, vol. 26, no. 4, p. 106399, Mar. 2023, doi: 10.1016/j.isci.2023.106399.
- [7] C. Zheng, Y. Sun, X. Ye, H. Chen, and H. Ji, “Establishment and characterization of primary lung cancer cell lines from Chinese population,” *Acta Pharmacol. Sin.*, vol. 32, no. 3, pp. 385–392, Mar. 2011, doi: 10.1038/aps.2010.214.
- [8] “Freezing Cells - US.” Accessed: Sep. 14, 2024. [Online]. Available: <https://www.thermofisher.com/us/en/home/references/gibco-cell-culture-basics/cell-culture-protocols/freezing-cells.html>
- [9] “Cell Culture | Parameters for Healthy Cells,” ibidi. Accessed: Sep. 14, 2024. [Online]. Available: <https://ibidi.com/content/435-parameters-for-healthy-cells>
- [10] “Introduction to CRISPR screening | IDT,” Integrated DNA Technologies. Accessed: Sep. 14, 2024. [Online]. Available: <https://www.idtdna.com/pages/education/decoded/article/overview-what-is-crispr-screening>
- [11] M. Selby *et al.*, “3D Models of the NCI60 Cell Lines for Screening Oncology Compounds,” *SLAS Discov.*, vol. 22, no. 5, pp. 473–483, Jun. 2017, doi: 10.1177/2472555217697434.
- [12] D. A. Cristaldi *et al.*, “A Reliable Flow-Based Method for the Accurate Measure of Mass Density, Size and Weight of Live 3D Tumor Spheroids,” *Micromachines*, vol. 11, no. 5, p. 465, Apr. 2020, doi: 10.3390/mi11050465.
- [13] Thermo Fisher Scientific, “generation-cancer-spheroids-tips-tricks-application-note.pdf,” Generation of cancer spheroids— tips and tricks. Accessed: Sep. 15, 2024. [Online]. Available: <https://assets.thermofisher.com/TFS-Assets/BID/Application-Notes/generation-cancer-spheroids-tips-tricks-application-note.pdf>
- [14] “Biosafety Levels and Their Meaning - Eppendorf US.” Accessed: Sep. 17, 2024.

- [Online]. Available:
<https://www.eppendorf.com/us-en/lab-academy/life-science/cell-biology/biosafety-levels-and-their-meaning/>
- [15] “1910.1030 - Bloodborne pathogens. | Occupational Safety and Health Administration.” Accessed: Sep. 17, 2024. [Online]. Available:
<https://www.osha.gov/laws-regs/regulations/standardnumber/1910/1910.1030>
- [16] “CFR - Code of Federal Regulations Title 21.” Accessed: Sep. 17, 2024. [Online]. Available:
<https://www.accessdata.fda.gov/scripts/cdrh/cfdocs/cfcfr/CFRSearch.cfm?fr=876.5885>
- [17] R. J. Geraghty *et al.*, “Guidelines for the use of cell lines in biomedical research,” *Br. J. Cancer*, vol. 111, no. 6, pp. 1021–1046, Sep. 2014, doi: 10.1038/bjc.2014.166.
- [18] J. L. Almeida, K. D. Cole, and A. L. Plant, “Standards for Cell Line Authentication and Beyond,” *PLoS Biol.*, vol. 14, no. 6, p. e1002476, Jun. 2016, doi: 10.1371/journal.pbio.1002476.
- [19] W. M. Bonner *et al.*, “ γ H2AX and cancer,” *Nat. Rev. Cancer*, vol. 8, no. 12, pp. 957–967, Dec. 2008, doi: 10.1038/nrc2523.
- [20] I. Padmalayam and M. J. Suto, “Chapter Twenty-Four - 3D Cell Cultures: Mimicking In Vivo Tissues for Improved Predictability in Drug Discovery,” in *Annual Reports in Medicinal Chemistry*, vol. 47, M. C. Desai, Ed., in Annual Reports in Medicinal Chemistry, vol. 47. , Academic Press, 2012, pp. 367–378. doi: 10.1016/B978-0-12-396492-2.00024-2.
- [21] N. Takahashi *et al.*, “3D Culture Models with CRISPR Screens Reveal Hyperactive NRF2 as a Prerequisite for Spheroid Formation via Regulation of Proliferation and Ferroptosis,” *Mol. Cell*, vol. 80, no. 5, pp. 828–844.e6, Dec. 2020, doi: 10.1016/j.molcel.2020.10.010.
- [22] J. R. Ferrarone *et al.*, “Genome-wide CRISPR screens in spheroid culture reveal that the tumor suppressor LKB1 inhibits growth via the PIKFYVE lipid kinase,” *Proc. Natl. Acad. Sci.*, vol. 121, no. 21, p. e2403685121, May 2024, doi: 10.1073/pnas.2403685121.
- [23] T. Stiff, S. Bayraktar, P. Dama, J. Stebbing, and L. Castellano, “CRISPR screens in 3D tumourspheres identified miR-4787-3p as a transcriptional start site miRNA essential for breast tumour-initiating cell growth,” *Commun. Biol.*, vol. 7, no. 1, pp. 1–11, Jul. 2024, doi: 10.1038/s42003-024-06555-1.
- [24] S. Gopal, A. L. Rodrigues, and J. S. Dordick, “Exploiting CRISPR Cas9 in Three-Dimensional Stem Cell Cultures to Model Disease,” *Front. Bioeng. Biotechnol.*, vol. 8, Jun. 2020, doi: 10.3389/fbioe.2020.00692.
- [25] C. Garcia-Canton, A. Anadón, and C. Meredith, “ γ H2AX as a novel endpoint to detect DNA damage: Applications for the assessment of the *in vitro* genotoxicity of cigarette smoke,” *Toxicol. In Vitro*, vol. 26, no. 7, pp. 1075–1086, Oct. 2012, doi: 10.1016/j.tiv.2012.06.006.
- [26] A. Ivashkevich, C. E. Redon, A. J. Nakamura, R. F. Martin, and O. A. Martin, “USE OF THE γ -H2AX ASSAY TO MONITOR DNA DAMAGE AND REPAIR IN TRANSLATIONAL CANCER RESEARCH,” *Cancer Lett.*, vol. 327, no. 1–2, pp. 123–133, Dec. 2012, doi: 10.1016/j.canlet.2011.12.025.

C. BME 400 Final Report

CRISPRi Screening in Cancer Spheroids to Investigate Factors in Genome Stability

December 17, 2024

Abstract

Current two-dimensional (2D) models have limited biological relevance to the tumor microenvironment, often failing to replicate the complex cell-cell interactions found *in vivo*. The objective of this project is to select a suitable cancer cell line to develop and refine a three-dimensional (3D) cell culture model that better recapitulates the tumor microenvironment. Toward this end, A549, a non-small lung cancer line (NSCLC) with mutations to tumor protein 53 and a relative sensitivity to the three main cancer drugs, was cultured and a doubling time of 23 hours was established. An optimal static spheroid formation protocol was engineered and tested by seeding spheroids with A549s at various densities. BioTek cytation images of these spheroids were taken and processed using ImageJ and MATLAB. From this data it was determined that ideal A549 aggregates were formed at a seeding density of 50k-75k.

In the future, the relevance of this information can be applied to scale up spheroid formation to prepare for γ H2AX—a sensitive marker for DNA double-strand breaks (DSBs)—staining and a high-throughput genome-wide CRISPRi screen. The team will provide valuable insights that bridges the gap between *in vitro* studies and clinical applications leading to more effective cancer therapies and improved patient outcomes.

1. Introduction

1.1. Motivation and Societal Impact

Cancer impacts millions of lives each year, with an estimated 2.2 million new cases and 736,790 deaths reported in the United States in 2024 [1], and not mentioning those still living with cancer. In the past, using 2D monolayer with CRISPR screening has provided extensive knowledge on what drives cancer cell high growth rate and what therapeutic treatments they are susceptible to. However, it fails to accurately represent the 3D tumor environment, leaving out important factors such as genome stability and DNA damage regulatory genes. Providing a replicated 3D tumor environment that is compatible with CRISPR screening can help identify sources of DNA damage not found with 2D cultures. This can lead to early detection, new therapeutic tools, and possibly genomic therapies.

1.2. Problem Statement

Although previous CRISPR screening in 2D monolayers has provided useful knowledge on cancer drivers and therapeutic susceptibilities, it lacks an element of biological relevance to an *in vivo* environment. Therefore, our team was tasked with developing a cell culture method that is compatible with a 3D environment and CRISPRi screening in order to identify sources of DNA mutations in the tumor environment. To achieve this goal, the team must select a viable cell line for the screen, create and optimize a spheroid formation protocol, and develop a protocol to stain for γ H2AX: a histone variant that is a sensitive marker for DNA damage.

2. Background

2.1. Background Research

2.1.a. 3D spheroids

The study of cancer and its progression has evolved significantly, with advancements in genetic and molecular biology providing deeper insights into the mechanisms underlying tumorigenesis. In this context, 3D spheroids have emerged as a superior model compared to

traditional 2D monolayer cultures. 2D cultures are limited in their ability to replicate the complex cellular interactions and microenvironment characteristic of actual tumors, which can lead to poor results in drug response and tumor behavior studies [2]. In contrast, 3D spheroids, composed of cancer cells aggregated in a spherical configuration, better mimic the architecture and growth phenotypes observed in tumors as seen in Figure 1 [2]. This structural arrangement facilitates critical cell-cell and cell-matrix interactions, which are essential for maintaining cellular functions and responses to therapeutic agents [3]. From these interactions, more DNA damage is observed in 3D spheroid cultures compared to 2D monolayers due to the hypoxic and nutrient-deprived microenvironment, which generates reactive oxygen species (ROS) and impairs DNA repair pathways [4]. Consequently, the use of 3D spheroids enhances the fidelity of experimental models, allowing for a more accurate investigation of tumor structure, cancer progression, and the efficacy of potential treatments.

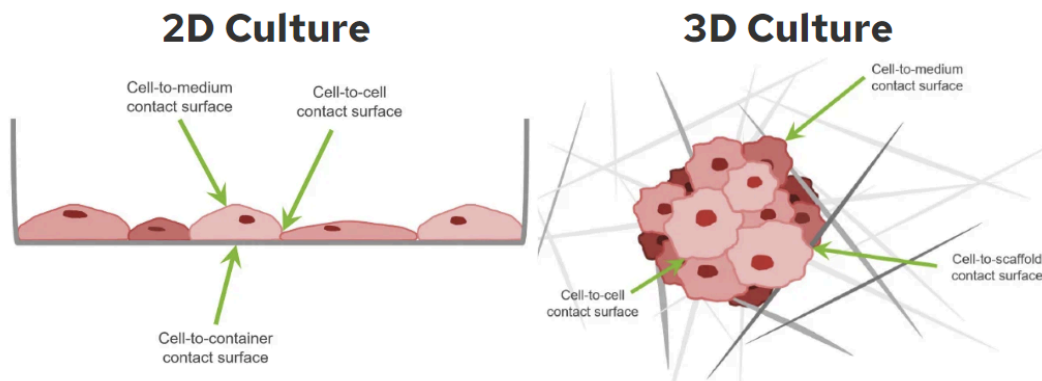


Figure 1: Comparison of 2D versus 3D tissue culture [5].

2.1.b. Cell line

Creating a cancer cell line is a multifaceted process that begins with the collection of tumor tissue from a cancer patient, typically obtained through biopsy or surgical resection. Once the tissue is collected, it undergoes processing to isolate individual cancer cells, often involving enzymatic digestion or mechanical disruption to break down the extracellular matrix and liberate the cells [6]. These isolated cells are then placed in a nutrient-rich culture medium designed to support their growth, containing essential nutrients, growth factors, and an appropriate pH. As the cells proliferate, certain cancer cells may demonstrate faster growth rates; these cells can be selected for further cultivation. Once cultivated, a cell line is then defined as continuous (infinite) or finite. For the purposes of this design project, both infinite and finite cell lines have been considered. After establishing a stable culture, the cell line undergoes thorough characterization to confirm its cancerous properties, assessing factors such as cell morphology, growth rate, and genetic markers. Finally, once a stable and well-characterized cell line is established, it is cryopreserved in liquid nitrogen for long-term storage [7]. The cell lines considered must be able to display a 3D tumor environment upon spheroid formation, be compatible with CRISPRi screening, and express γ H2AX.

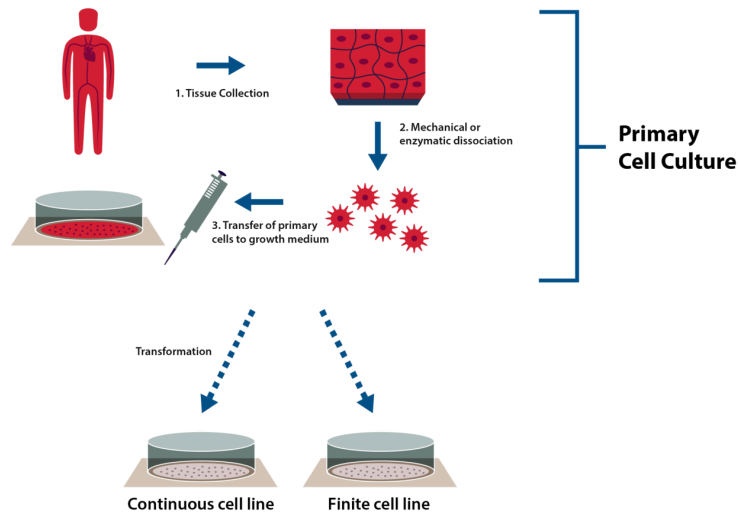


Figure 2: Cell line creation [6].

2.1.c. CRISPRi Screening

CRISPR interference (CRISPRi) is a powerful tool for gene regulation that allows for the targeted repression of genes without directly inducing double-strand breaks in the DNA. This technique utilizes a catalytically inactive Cas9 (dCas9) protein fused to transcriptional repressors, enabling precise control over gene activity [8]. In this project, lentiviral vectors are employed to deliver the CRISPRi components into the target human lung cancer cell line, which is A549. The use of lentivirus is advantageous due to its ability to efficiently transduce a wide range of dividing and non-dividing cells, ensuring stable expression of the dCas9 and guide RNA (sgRNA) constructs necessary for effective gene repression [8]. For instance, SOX2 is a transcription factor associated with stemness and tumor progression in non-small cell lung cancer (NSCLC). It regulates self-renewal, differentiation, and pluripotency, making it a key driver of cancer stem cell (CSC) properties. SOX2 is frequently co-expressed with surface markers such as CD133 and CD44, which are commonly used to identify CSCs. CD133, a marker of tumorigenicity, is often upregulated in SOX2-expressing cells and is linked to enhanced self-renewal and tumor initiation [9]. Similarly, CD44, which plays a role in cell adhesion and migration, is associated with SOX2 in promoting metastasis and resistance to therapies [9]. These surface markers, in conjunction with SOX2, delineate a subpopulation of CSCs that contribute to lung cancer progression, recurrence, and therapeutic resistance, making them a target for gene therapies in NSCLC.

By applying CRISPRi in a 3D spheroid model, the team's design aims to better mimic the in vivo environment for the screen, which is essential for understanding the complex interactions that drive cancer progression and treatment resistance. As stated, this approach facilitates the identification of gene targets that contribute to DNA damage responses and cellular survival in a more physiologically relevant context compared to traditional 2D cultures.

2.1.d. γ H2AX

γ H2AX is a phosphorylated variant of the histone H2AX, which serves as a sensitive marker for DNA double-strand breaks (DSBs) [2]. The phosphorylation of H2AX occurs rapidly in response to DNA damage, making it an essential marker for monitoring the integrity of the

genome [3]. In cancer research, γ H2AX staining is utilized to assess the extent of DNA damage and the efficacy of DNA repair mechanisms.

γ H2AX allows for the quantification of DNA damage induced by various genetic suppression within the 3D spheroid model. γ H2AX staining is more challenging in a 3D environment due to the compact and multilayered structure of spheroids creating physical barriers that hinder the diffusion of antibodies into the inner core. This lack of uniform antibody penetration can result in incomplete staining, where only the outer layers of the spheroid are adequately labeled, leading to a skewed representation of DNA damage throughout the entire structure.

If the γ H2AX staining is done correctly in a CRISPRi screen within 3D spheroids, it enables accurate assessment of DNA damage across the entire structure, including inner cell layers, providing a comprehensive view of the effects of gene knockdowns on genome stability. This ensures that the screening results are not biased toward cells in the outer layers and allows for the identification of key genetic contributors to DNA repair and damage pathways in 3D. The response of dsDNA breaks can effectively evaluate how specific gene targets contribute to genomic stability, tumor biology, cancer progression, and therapeutic vulnerabilities.

It is important to choose a cell line that has been proven to be sensitive (or has a baseline of sensitivity) to DNA damage for successful γ H2AX staining. One way to determine this criteria is to see whether that cell line is sensitive to drugs that can induce DNA damage. The drugs chosen for this purpose include cisplatin, oxaliplatin, and bleomycin.

Cisplatin is a platinum-based chemotherapeutic agent widely used in the treatment of various solid tumors, including lung, ovarian, and bladder cancers. Its mechanism involves the formation of DNA cross-links, which interfere with DNA replication and transcription, leading to cell cycle arrest and apoptosis [10]. Another platinum-based compound with similar anticancer mechanisms is oxaliplatin, which is more effective in some forms of cancer such as colorectal cancer [11]. On the other hand, bleomycin, an antitumor antibiotic, induces DNA damage through the generation of free radicals, leading to both single and double-strand DNA breaks. Unlike platinum-based drugs, bleomycin's mechanism of action makes it suitable for combination therapy, targeting cancer cells through a different pathway [12].

Drug sensitivity is quantitatively measured by Z-score, which measures drug sensitivity of a specific cell line to a specific drug in comparison to other cell lines. Negative Z-scores reflect sensitivity, with lower values indicating greater susceptibility to the drug's cytotoxic effects, while positive Z-scores indicate resistance [13]. A Z-Score of ± 1 is considered significant [14]. These differential responses across cancer types and drug classes highlight the complexity of cancer treatment and the need for tailored therapeutic approaches based on specific tumor characteristics.

2.1.e. *Previous Literature*

In previous research, 3D cancer spheroid models have been utilized for CRISPR-Cas9 screens with and without using an exogenous extracellular matrix (ECM).

Takahashi *et al*, 2020 utilized 4% Matrigel, an ECM based approach, as bio-scaffolding to plate A549 cell spheroids. They concluded that their seeded spheroid models *in-vivo* interactions between cells and the ECM were more biologically relevant than former 2D studies [15]. The study's CRISPRi library gene list as well as its phenotypes have a strong overlap with those from Han *et al*, 2020, suggesting that inner cells of ECM-based 3D spheroids may be similar to ECM-free 3D spheroids.

Han *et al*, 2020 devised an ECM-free scalable method to propagate 3D spheroids [16]. H23 cells were seeded on pre-treated ultra-low attachment or poly-HEMA-covered plates in RPMI 1640 medium with 0.75% methylcellulose. Similar methods have been found to be implemented in other studies, such as by Ferrarone, *et al* (2024) and by Stiff, *et al* (2024) [17], [18]. A total of around 210,000 single guide RNAs (sgRNAs) were designed for the CRISPRi library with around 10 sgRNAs per gene [16]. Samples were used to calculate growth and/or tested for drug resistance [16].

3D stem cell models have also been genomically edited via CRISPR-Cas9 to better model various disease phenotypes. Cells can be assembled with or without supporting scaffolds. Some example scaffolds include laminin-rich Matrigel, collagen matrix, or those that are naturally derived from decellularized tissues [19].

2.2. Client Information

The clients, Ms. Carley Schwartz and Dr. Gaelen Hess are researchers at the Hess Laboratory within the Wisconsin Institutes for Medical Research (WIMR). The group's research focus is on high-throughput functional genomics to investigate DNA repair and pathogenic effectors [20]. Dr. Hess is an assistant professor in the Department of Biomolecular Chemistry and Center for Human Genomics and Precision Medicine [20]. Ms. Carley Schwartz is the lab manager and a research intern at the Hess Laboratory.

2.3. Design Specifications

The project aims to develop a 3D cell culture method compatible with CRISPRi screening to identify sources of DNA mutations in tumors, utilizing γ H2AX as a marker for DNA damage. The client requires the team to select an appropriate human lung cancer cell line, to which A549 and NCI-H23 meet that criteria [21]. This protocol must focus on optimal seeding density, targeted at 50,000 cells/cm², and viscosity adjustments through methylcellulose to ensure uniform spheroid size and maximal cell growth [16]. Additionally, a γ H2AX staining protocol needs to be established following the guidelines of scalability, and be reproducible with optimized antibody concentrations and reagent volumes.

With the use of human contaminants, safety measures are to adhere to all Biosafety Level 2 (BSL-2) standards, implementing containment protocols and proper waste disposal [22]. The accuracy and reliability of the spheroid formation process will be monitored through live-cell imaging to confirm a ~30% cell death rate, and at least three staining experiments will be conducted to establish statistical confidence in the results.

For life in service, cells will require approximately two months in culture for protocol optimization, spheroid formation following the culture, and a minimum of six days post-infection with the CRISPRi library. To ensure long-term viability, cells will be cryopreserved at temperatures below -135°C. Operating conditions will be maintained at 37°C, 5% CO₂, and 90–95% humidity [23].

The design will prioritize ergonomic efficiency in the experimental setup, aiming to enhance throughput while accurately mimicking *in vivo* conditions. Spheroids are targeted to be ~500 μ m in diameter, and accurate assessments of mass and density will be facilitated through specialized equipment [24]. The project budget is capped at \$1,000, with a focus on utilizing existing lab materials to keep costs manageable.

Lastly, the team will conduct a competitive analysis of existing 3D cancer models, focusing on methodologies that incorporate ECM-free approaches. Key milestones include cell

line creation and protocol development by October 2024, optimization of spheroid formation by November 2024, completion of the γ H2AX staining protocol by December 2024, execution of high-throughput CRISPR screening by February 2025, followed by analysis of phenotypic differences before May 2025. For more detailed information see Appendix 12.1.

3. Preliminary Designs

3.1. Cell Lines

3.1.a. NCI-H23

NCI-H23 is an epithelial-like cell that was isolated from the lung of a Black, 51-year-old, male patient with adenocarcinoma. H23 is a non-small cell lung cancer (NSCLC) with a doubling time 38 hours [25]. Although H23 is considered an adherent cell line, it is known to have weaker cell-cell interactions due to the lack of E-cadherin at the cell membrane [26]. H23 exhibits TP53, KRAS-12, epidermal growth factor receptor (EGFR), L-dopa decarboxylase-negative, and a high degree of c-myc DNA amplification (20-fold) but no detectable amplification of c-myc RNA [25]. Additionally, these cells have a reported colony forming efficiency of 9.7% in soft agarose [2]. H23 is very sensitive to bleomycin (Z-score average of -1.22), somewhat sensitive to cisplatin (-0.73), and somewhat resistant to oxaliplatin (+0.41) [27]. One vial of these cells from ATCC is \$555.00 [25].

3.1.b. A549

A549 is an adherent cell line isolated from the lung tissue of a White, 58-year-old male with lung cancer [28]. This NSCLC adenocarcinoma has a doubling time of 22 hours [28]. A549 exhibits mutations for TP53, EGFR, PIK3CA, ALK, and PTEN [29]. A549 is very sensitive to bleomycin (Z-score average of -1.04), somewhat sensitive to cisplatin (-0.12), and very sensitive to oxaliplatin (-1.45). One vial of these cells from ATCC is \$555.00 [28].

3.1.c. OVCAR-5

OVCAR-5 is a cell line established from ascites fluid from a 67 year old, female, non-treated patient with a metastatic gastrointestinal cancer that has presented as advanced-stage ovarian carcinoma [30]. The doubling time for this cell line is 27 hours [31]. OVCAR-5 expresses mutations for KRAS-12, EGFR, and CLAUDIN-4 [30]. This ovarian cancer cell line is somewhat resistant to bleomycin (Z-score average of +0.67), somewhat resistant to cisplatin (+0.52), and somewhat sensitive to oxaliplatin (-0.18) [32]. This cell line is not currently available for the team to purchase [30], [33].

3.2. Spheroid Formation Protocol

3.2.a. Treated Tissue Culture Plates

Treated tissue cultures method is a scaffold-free method [34]. It uses low attachment plates, which are coated with a hydrophilic polymer (such as hydrogel or poly-HEMA) to prevent specific and nonspecific cell adhesion to the culture vessel [34]. Thus, high cell-to-cell interactions allow cells to cluster together, form cell aggregates and eventually into 3D spheroids (Figure 3) [35]. The number and size of spheroids formed depend on the shape of the wells, with U-bottom wells creating a single and well-defined spheroid for each well (Figure 3) [35]. The whole process usually takes up to 4 days [35].

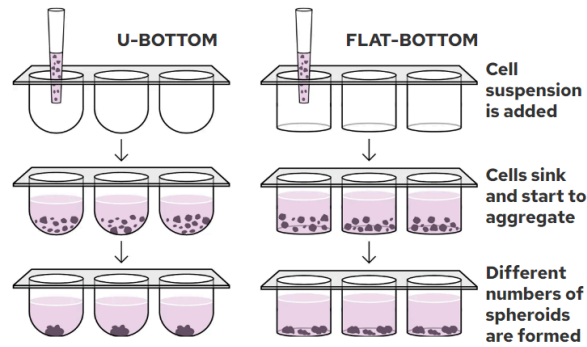


Figure 3: Spheroid formation by treated tissue culture plates/low attachment plates with U-bottom wells (left) and flat-bottom (right) [34].

3.2.b. Hanging Drop

Similar to the treated tissue cultures method, the hanging drop method is also anchorage-independent [34]. Currently, two methods to form hanging drops exist. The first method is the traditional method, involving dispensing small amounts of cell suspension on a surface, usually the underside of a petri dish, then flipping the dish upside down to form hanging drops (Figure 4) [36]. However, this method has a high risk of droplet fusion (if space is limited), loss of droplet shape (because of frequent media change), and droplet evaporation [37], [38]. The second method using hanging drop plates negates these risks [38]. Each plate consists of a lid to maintain sterility and a main component, which has a water reservoir to prevent evaporation and access holes on top where cells can be added to form hanging drops on the bottom (Figure 5a) [38]. Hanging drops create a flexible spheroid cell-to-liquid interface, allowing cells to gravitationally aggregate and form because of high cell-to-cell interactions (Figure 5b) [34]. The whole process usually takes up to 4 days [35].

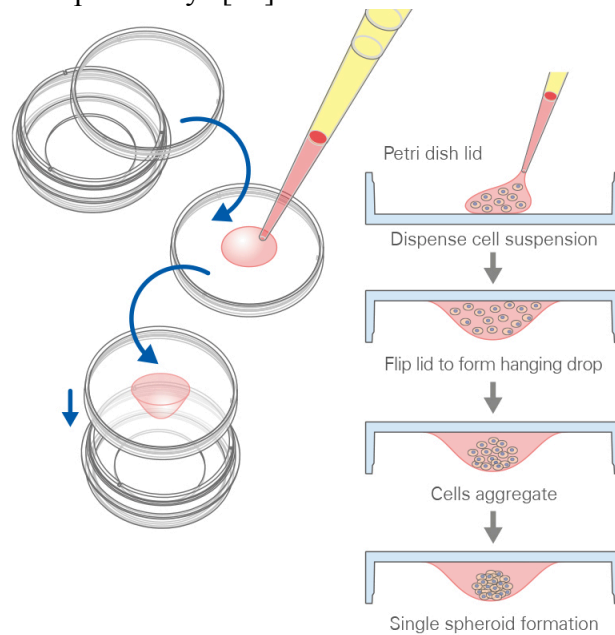


Figure 4: Spheroid formation by flipped hanging drop method [36].

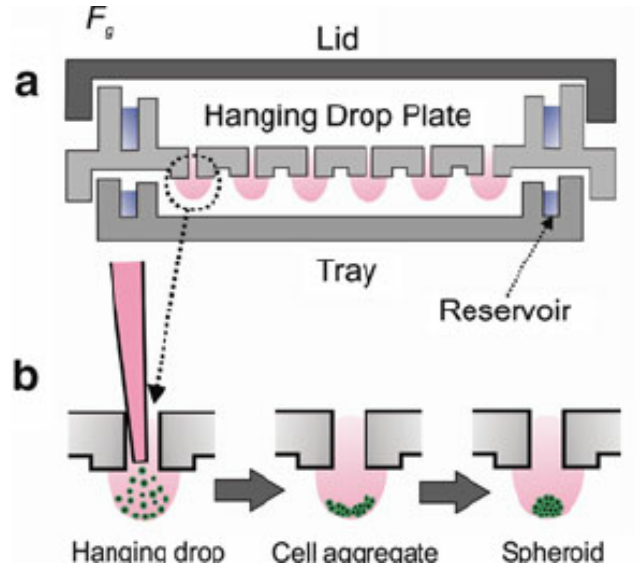


Figure 5: Spheroid formation by hanging drop plate [34].
 a. Composition of a hanging drop plate, b. Spheroid formation.

3.2.c. Matrigel

The Matrigel method is a scaffold-dependent method [34]. Cells are cultured in low attachment plates for 4 days to form aggregates, then Matrigel is added to encourage cell-to-ECM interactions which will form spheroids after another 2 days (Figure 6) [39].

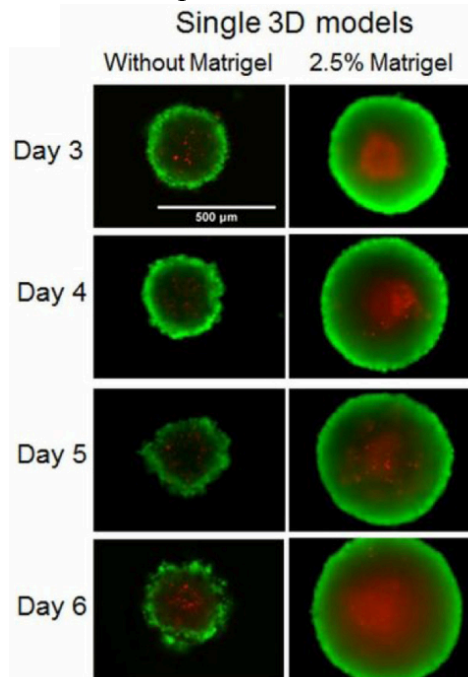


Figure 6: Comparison of single 3D spheroids grown without and with Matrigel [40].
 Green indicates live cells and red indicates dead cells.

4. Preliminary Design Evaluation

4.1. Design Matrix for Cell Line

Table 1. Design matrix for Cell Line.

Design Criteria (weight)	NCI-H23		A549		OVCAR-5	
	Score X/5	Weighted score	Score X/5	Weighted score	Score X/5	Weighted score
Viability (20)	4	16	5	20	3	12
Adhesion (20)	3	12	4	16	4	16
Reproduction Speed (20)	3	12	5	20	4	16
Drug Sensitivity (15)	4	12	5	15	3	9
Genetic Mutation (15)	5	15	3	9	2	6
Ease of Procurement (10)	5	10	5	10	5	10
Total Score (100)	77		90		69	

Determination of Criteria and Weights & Justification of Assigned Scores: After discussion about the important criteria to evaluate the initial designs against our Product Design Specifications, the following categories were chosen.

Viability (20):

Cell viability is a measure of the proportion of live, healthy cells within a population [41]. For this category, viability is a measure of how robust the cell is and how well it can withstand various stresses. The highest ranked cell line is predicted to be the healthiest after fluid shear stress, transduction, and cisplatin treatment.

NCI-H23 cells were assigned a score of 4/5 for its and >20% cell death after 72 hours when subjected to [3 μ M] of cisplatin and its sensitivity to cisplatin treatment [42]. The NCI H23 cell line showed a higher expression level of DNA repair proteins after cisplatin treatment compared to A549 cells. In previous experiments for both NCI-H23 and A549, a negative correlation between cell viability and DNA damage induction upon cisplatin treatment has been noted [42]. A549 cells were assigned a score of 5/5 for its >20% cell death after 72 hours when subjected to [3 μ M] of cisplatin and its greater resistance to cisplatin treatment compared to NCI-H23 cells [42]. OVCAR-5 scored a 3/5 for its severe morphological changes (decreased spheroid-forming capacity, reduced cytoskeleton organization, and chromosomal instability) after exposure to <1 dyne/cm² fluid shear stress and >50% cell death after 72 hours when subjected to [2.5 μ M] cisplatin [43], [44].

The type of transduction done to these three cell types may either help or hinder the cell health depending on the transduction virus used (lentivirus or siRNA) and what it is used in

conjunction with (Nrf2, cisplatin, or both). The trend in scientific literature indicates that these three cell types have the same viability, within $7\pm 3\%$, when transfected under the same conditions [42], [43], [44]. Since the vectors encoding the components necessary for CRISPR/Cas are large, they result in low transduction efficiency and cell viability [45]. To overcome those obstacles, the team can add exogenous small plasmids that increase transduction efficiency up to 40-fold and cell viability up to 6-fold [45].

Adhesion (20):

This criteria refers to the ability of cells to attach to the surface of the culture vessel or substrate they are grown in, as well as to attach to other cells. For this criteria, cell lines were therefore judged on whether they were considered an adherent cell line and on their strength of cell-cell interactions. Because adhesion is an important factor in cells being able to effectively form spheroids, this criteria was given a high weight for scoring.

Although NCI-H23 is considered an adherent cell line [25], it is known to have weaker cell-cell interactions due to the lack of E-cadherin at the cell membrane [46]. Additionally, NCI-H23 has presented difficulty in forming spheroids in previous literature, particularly using the Matrigel method [26]. For these reasons, NCI-H23 scored 3/5. Similarly, although A549 and OVCAR-5 are considered adherent cell lines [28], [47], [48], it is also known to have a lack of E-cadherin and thus weaker cell-cell interactions [46], [49]; thus, A549 and OVCAR-5 scored 4/5.

Reproduction Speed (20):

Reproduction Speed is the time it takes the cells to grow and reproduce. In order to look specifically at the cell line, and exclude any factors caused by the spheroid formation protocol, the team used the cell lines' doubling speed. The doubling speed is the amount of time it takes for a cell population to double in number. Essentially, it is how quickly a cell line replicates itself in culture [50].

As each cell line has a range of doubling times depending on the culture medium used, the most typical doubling times were used to determine this category's score. With the longest doubling time of about 38 hours [25], [51], NCI-H23 scored the lowest at a 3/5. A549 is the fastest of the three cell lines with a typical doubling time of about 22 hours [28], [51], scoring a 4/5. Lastly, OVCAR5 has a doubling time of 27 hours [31], [51] scoring a 4/5.

Drug Sensitivity (15):

Gamma-H2AX foci staining is used for this project to assess the extent of CRISPR-Cas9-based DNA damage. Cell lines that are particularly sensitive to drugs that induce DNA damage will have a baseline of DNA damage to refer to, making them appropriate for this project.

Z-scores for NCI-H23, A549, and OVCAR-5 are obtained from the Genomics of Drug Sensitivity in Cancer Project [52]. NCI-H23 has a Z-score average of -1.22 for bleomycin, -0.73 for cisplatin and +0.41 for oxaliplatin [27]. A549 has a Z-score average of -1.04 for bleomycin, -0.12 for cisplatin, and -1.45 for oxaliplatin [53]. OVCAR-5 has a Z-score average of +0.67 for bleomycin, +0.52 for cisplatin, and -0.18 for oxaliplatin [32]. A549 has 2 significant negative Z-scores, scoring a 5/5. NCI-H23 has a negative Z-score for 1 of the drugs, scoring a 4/5. OVCAR-5 has none, scoring a 3/5.

Genetic Mutation (15):

Analysis of genetic mutations inherent to various cell lines is crucial in determining the ideal cell line for our project. A high score in this category means genetic mutations that are helpful must work with for any step in the adhesion process, Cas9 attachment, CRISPR edit, or screen. Using Han's ten selected lung cancer cell lines, a set of ideal cell line mutations were created including: TP53 (tumor protein p53), a gene that encodes a tumor suppressor protein containing transcriptional activation, DNA binding, and oligomerization domains [54]; EGFR, a cell surface protein that binds to epidermal growth factor leading to cell proliferation [55]; and KRAS, Kirsten rat sarcoma viral oncogene homolog causes malignancies based on an error encoding a protein that is a member of the small GTPase superfamily [56]. These mutations may either cause or enhance the malignancy of existing cancer cells. Since TP53 is the most frequently mutated gene in human cancer, the presence of this mutation is crucial for biological relevance to half of all cancer cases [57].

NCI-H23 has all of the key mutations including TP53, KRAS, EGFR, a high degree of c-myc DNA amplification, and no counterproductive mutations, so it was given the highest score of 5/5 [25]. A549 scored second highest with 3/5 because it contains EGFR and TP53 which are key, but also PIK3CA, ALK, and PTEN which aren't key but won't impede the success of the project [29]. The protein encoded by PIK3CA, phosphatidylinositol 3-kinase, is oncogenic and represents the catalytic subunit, which uses ATP to phosphorylate PtdIns, PtdIns4P and PtdIns(4,5)P2 [58]. ALK, anaplastic lymphoma receptor tyrosine kinase, a chromosomal rearrangement which results in creation of multiple fusion genes in tumorigenesis. PTEN guards the genome by controlling multiple processes of chromosome inheritance [59]. Although OVCAR-5 has KRAS, a key mutation, it was ranked the lowest at 2/5 due to the lack of essential protein TP53 and the presence of additional mutation CLDN4, which is a mutation of integral membrane proteins [30], [60], [61].

Ease of Procurement (10):

This criteria assesses how easy it is to obtain the cell line. It is expensive to purchase an entirely new cell line from ATCC, Sigma-Aldrich or similar sources. Therefore, the chosen cell line should be either already available in Hess Lab or can be obtained from other labs, preferably within Wisconsin Institutes for Medical Research (WIMR); however, that process can be time-consuming.

All 3 chosen cell lines are readily available in Hess Lab, both wild-type and genetically-modified, and therefore all scored a 5/5.

4.2. Design Matrix for Spheroid Formation

Table 2. Design matrix for Spheroid formation.

Design Criteria (weight)	Treated Tissue Culture Plates		Hanging Drop		Matrigel	
	Score X/5	Weighted score	Score X/5	Weighted score	Score X/5	Weighted score
Ease of Fabrication (30)	5	30	4	24	3	18
Scalability (20)	4	16	3	12	5	20
Uniform properties (20)	4	16	3	12	2	8
Cost (15)	3	9	4	12	2	6
Safety (15)	5	15	5	15	5	15
Total Score (100)	86		75		67	

Determination of Criteria and Weights & Justification of Assigned Scores: After discussion about the important criteria to evaluate the initial designs against our Product Design Specifications, the following categories were chosen.

Ease of Fabrication (30):

This criteria assesses the complexity of the 3D spheroid formation protocol, which is one of the important factors when choosing an appropriate fit-for-purpose 3D culture model [62]. Specifically, this criteria assesses whether additional specialized training or equipment will be needed for the protocol, on top of those necessary for 2D traditional cell culture methods. Also, the time needed to complete the protocol once will also be considered, as it can directly affect the length of the project due to the large number of cells needed for genome-wide CRISPR screening.

Treated tissue culture plates method scores a 5/5 as the cell seeding method is not particularly labor-intensive, only requires 1 media change and takes the least time (up to 4 days). Hanging drop method also takes the same amount of time but the cell seeding method requires more careful technique to ensure hanging drop viability and daily media change, thus scoring a 4/5. The Matrigel method does not require extensive cell seeding technique and daily media change but requires longer time (6 days), scoring a 3/5.

Scalability (20):

This criteria assesses the plausibility to use these methods to generate a large amount of viable spheroids. A high-throughput method is preferred because, per the client, each CRISPR screen would require 100-150 million cells per replicate to generate a sufficient amount of spheroids [16]. Because a high-throughput method is very important to successfully completing a spheroid CRISPR screen, this criteria was given a high weight for scoring.

The treated tissue culture plates method is known for its high-throughput capabilities because treated plates have a non-adherent coating that promotes rapid cell aggregation into spheroids [63]. This method also allows for a greater amount of spheroids to be made in a single plate because it can accommodate large volumes of media [64]. For these reasons, this method was given a score of 4/5 for scalability. The Hanging Drop method, in contrast, is considered to be low-throughput in nature. One reason for this is because the method requires a lot of manual handling and is prone to error when setting up and transferring the drops [65]. Additionally, scalability is limited with the Hanging Drop method because each tissue culture plate can only hold a limited amount of drops, and each drop only tends to form one spheroid [37]. For these reasons, this method was given a score of 3/5 for scalability. Lastly, the Matrigel method is considered high-throughput because, like the treated tissue culture plates method, it can allow for a comparatively larger volume of media to generate spheroids compared to the Hanging Drop method. Additionally, this method involves polymerizing a layer of Matrigel matrix onto tissue culture plates to promote spheroid formation [26], and unlike the other two methods, has been shown to improve the growing environment of spheroids due to its mimicry of in-vivo tumor conditions [40]. For these reasons, this method was given a score of 5/5 for scalability.

Uniform Properties (20):

This criteria outlines the reliability of the spheroid formation to produce uniform properties. The ideal formation protocol has no batch-to-batch variability. Minimal variation in size, shape, and porosity is expected so the experiment will have a higher chance of reproducibility and fewer outliers. It is important to note that research has shown that differences in spheroid formation arise from the differences in cell-biomaterial interactions rather than due to differences in cell viability or proliferation [66].

The treated tissue culture plates scored the highest in this category, 4/5, due to their minimal variation between batches. Additionally, treated tissue culture plates have a larger amount of spheroids per batch (see *Scalability* section), making it easier to harvest a uniform size or shape [62]. Due to the low-throughput nature of Hanging Drop, making it harder to harvest uniform spheroids, and due to its 10% to 15% size variation between drops, this method scored a 3/5 [67]. Matrigel, animal-based basement membrane extract, has a batch-to-batch variation that can hinder reproducibility, so a score of 2/5 was assigned to it [66]. Cells seeded in matrigel also tend to form large cell aggregates with looser morphology rather than spheroids with round, regular shape [66].

Cost (15):

This category is scored based on the expenses of the materials needed to perform the spheroid formation protocol. While the client stressed the importance of a working protocol they also expressed their desire to keep costs to a minimum, as to be able to repeat the protocol and perform it on a larger scale.

Explain scoring:

Treated tissue culture plates scored a 3/5 in this category. Generally, the culture plates range from \$0.22-3.78 per well [68], with a need for an abundance of wells. The cost of treated tissue culture plates can vary depending on the manufacturer, the number of wells needed, and the volume of cell culture media required. The Hanging Drop method is rather inexpensive as it does not require specialized equipment or reagents. At around \$4 [69], [70] per 100 mm plate the Hanging Drop method is the most cost effective method out of all the options and scored a 4/5

[71]. Lastly, Matrigel scored a 2/5 in this category as the Matrigel itself is rather expensive, at \$380 for 10 mL of Matrigel [66], [72].

Safety (15):

This criteria assesses the ability of the method to be used in a laboratory environment without posing significant health risks to the user or the experimental results. This includes evaluating the potential hazards associated with the materials/reagents and processes used, the ease of handling, and the effectiveness of protective measures (e.g., PPE, ventilation). A method should be dependable, minimizing the risk of contamination or adverse effects that could compromise the safety of the personnel and the integrity of the experiments.

All of the methods proposed were able to obtain a 5. This is due to none of the methods using harmful reagents in the process of fabrication. All of the cell lines pose similar contamination risk but can be lowered considerably with proper PPE use.

4.3. Proposed Final Design

A549 was the cell line selected based on the team's decision matrix (Table 1). A549 is an adherent cell line which exhibits the most useful mutations (EGFR, TP53, PIK3CA, ALK, and PTEN) [28], [29]. Since this cell line has a doubling time of approximately 22 hours, the team will need to passage the cells 2-3 times a week [73]. Hess lab recommends a maximum of 20 passages be done on each cell line, so, once thawed, each cell line will have approximately 7 weeks before it must be bleached and another cell line thawed [73].

Treated tissue culture plates, or low attachment plates method was the selected spheroid formation protocol selected based on the team's decision matrix (Table 2). The process includes using hydrophilic polymer-coated plates [34], which are not labor intensive to make, and spheroid formation takes around 4 days [35], which is not time-consuming. Furthermore, as this method does not require a scaffold, it does not run the risk of high batch-to-batch variability due to using animal-derived matrices [34].

5. Design Process

5.1. Doubling time

5.1.a. Calculations used

See Appendices 12.2 and 12.3 for the full cell line maintenance and passaging protocols. Doubling time of cultured A549 cells was determined based on a Monday-Wednesday-Friday passaging schedule using equation 1 below:

$$t_d = \frac{\ln(2)}{\mu} \cdot 24h$$

Equation 1. Equation for calculating a cell line's doubling time. t_d = doubling time (days) and μ = growth rate (1/day). μ was obtained from the CytoFLEX cell confluency reading [74].

Once doubling time was determined, the volume of cells to passage was calculated using equation 2 below:

$$\frac{\# \text{ cells in culture solution}}{10 \text{ mL}} = \frac{\# \text{ cells needed to passage}}{\text{volume needed to passage}}$$

Equation 2. Equation for calculating the required volume for passaging. The number of cells needed to passage was estimated by looking at previous passage values and their resulting doubling times.

5.1.b. CytoFLEX Protocol

Comprehensive CytoFLEX protocol can be found below and in Appendix 12.3. Briefly, ensure daily clean has been run, load a well vortexed sample, alter settings to fast, click run, auto set axes from properties menu, gate a representative live cell population, and save file to BME team folder.

5.1.c. Final established doubling time

The cultured A549 cells' average established doubling time (for cells derived from second cryovial thawed) was 43.94 hours. When cutting down to 1.75 million cells for a two-day passaging period (Monday to Wednesday; Wednesday to Friday), the doubling time was closer to the doubling time of 22 hours observed in literature [28], averaging at 22.93 hours. However, when cutting down to below 1.75 million cells for a 3-day passaging period (Friday to Monday) the doubling time was much higher, averaging at 52.95 hours. It is also important to note that there was significant variability in measured doubling times with cells derived from the first cryovial of the semester, with some doubling times reaching up to 72 hours. A suspected reason for this variability was that step 8) in Appendix 12.3 was not followed consistently by all team members. This likely caused the CytoFLEX cell confluency reading to be underestimated for these passages and, as a result, the doubling times to be higher than expected. See Appendix 12.4 for the cell line maintenance benchling data.

5.2. Spheroid Protocol Selection

Once treated tissue culture plates were selected for our fabrication method, additional research was done to compare the client's spheroid formation protocol to protocols from other literature (section 2.1e). With the client's assistance, the selection of techniques was narrowed down to K. Han *et al.* and J. R. Ferrarone *et al.*. Both protocols were fairly similar in terms of complexity and required materials, but differed greatly in terms of time and scalability (table 3). Overall, a modified version of K. Han *et al.* was used for preliminary spheroid formation and seeding density testing. The main modifications made include using DMEM + 10% Fetal Bovine Serum (FBS) and 1% penicillin-streptomycin rather than the RPMI, using different TC plates, different poly-HEMA and Methylcellulose preparation methods, and exact seeding density. K. Han *et al.* used RPMI to prevent excessive aggregation of cells in the spheroid, but that should not be a concern for the team in this early testing phase. Other alterations made to the protocols will be discussed in the testing and discussion sections respectively. Deciding on an exact preliminary protocol to follow was essential for materials selection and ordering.

Table 3. Additional design matrix for spheroid formation protocols.

Design Criteria	1: K. Han <i>et al.</i>	2: J. R. Ferrarone <i>et al.</i>
Cost	Medium	Low
Spheroid Formation	Overnight*	24 hours

Time		
Scalability -cells/spheroid -spheroid #	20k-150k cells/cm2	1k cells/well
Required machinery & materials	1) 500uL media containing 0.75% methylcellulose 2)RPMI 1640 growth medium 3)96-well poly-HEMA TC plates	1) 80 µL of methylcellulose-containing medium 4)RPMI 5) V-bottom 96-well plates (low-attachment)

*time does not include Methylcellulose, poly-HEMA, media, or plate preparation.

6. Methods

6.1. Materials

Materials listed are required for culturing and passaging A549 cells and are already available in Hess Lab and thus require no purchasing. One 1 mL cryovial of A549 can be used up to 20 passages and 3 passages will be performed for each week of the project. In total, two cryovials of A549 were used, and a total of 20 passages (15 for the first cryovial and 5 for the second) were performed.

Table 4. Materials List.

Material	Part Specification	Quantity	Price
Already available in Hess Lab, requires no purchasing			
T-75 Flasks	Catalog #156800 - Thermo Fisher Scientific	1 for up to 20 passages	\$381.65 / case of 100 [75]
A549	CCL-185 - ATCC	N/A	\$555 / vial [28]
DMEM (high glucose)	Catalog #11965118 - Thermo Fisher Scientific	1x500 mL bottle 10 mL / passage (used in solution with Pen/Strep and FBS)	\$264.65 / 10 x 500 mL [76]
Penicillin- streptomycin (Pen/Strep)	Catalog #15070063 - Thermo Fisher Scientific	1% in solution with DMEM and FBS	\$23.65 / 100 mL of 5,000 U/mL [77]
Fetal Bovine Serum (FBS), Value	Catalog #A5256701 - Thermo Fisher Scientific	10% in solution with DMEM and Pen/Strep	\$314.37 / 500 mL [78]
PBS, pH 7.4	Catalog #10010072 - Thermo Fisher Scientific	2 mL / passage	\$244.00 / 6 x 1000 mL [79]

Trypsin-EDTA (0.05%), phenol red	Catalog #25300120 - Thermo Fisher Scientific	3 mL / passage	\$287.65 / 20 x 100 mL [69]
VWR® Flat-bottom 96-well polystyrene assay plates	Catalog #76446-962 - VWR	2 plates	\$494.24 / 100 plates [80]
Ethanol	Catalog #BP28184 - Fisher Scientific	33 mL	\$1,595.00/ (4 x 4L) [81]
Requires purchasing			
Poly-HEMA	Catalog P3932-10G - SigmaAldrich	1.3 g	\$238 / 10 g [82]* 1.3 g = \$30.94
Methylcellulose	Catalog M0512-100G - SigmaAldrich	2 g	\$51.40 / 100 g [83] * 2g = \$1.03
Accutase	Catalog A1110501 - ThermoFisher	Has not been used	\$60.65 / 100 mL [84]
Total:			\$92.62

6.2. Protocols

6.2.a. *Hess Lab Cell Line Maintenance*

Thawing Protocol: This protocol is developed by Hess Lab [85]. New media is prepared in a 15 mL tube. Slowly mix the media in the cryovial to thaw cell pellets. Once thawed, cells are pelleted by centrifugation and resuspended in fresh media. Cells need to be grown for typically 3 days before passaging further. See Appendix 12.2 for detailed protocol.

Passaging Protocol: This protocol is developed by Hess Lab [85]. This assumes passaging in a T75 flask (is using a T150, double the volume). Media is removed via aspiration. Wash cells with PBS then detach cells from the culture vessel with trypsin, then media is re-added to inactivate the trypsin. Cells are moved to a 15 mL conical tube, pelleted by centrifugation, and resuspended in new media. Add cell suspension to a new flask and add new media, enough to cover the cells. See Appendix 12.2 for detailed protocol.

6.2.b. *Spheroid formation protocol*

Spheroid formation protocol is based on Han *et al*, 2020 [16]. This semester, cells will be seeded at multiple densities to determine the optimal density for high spheroid numbers, uniform spheroid size, and high percent viability (see section 7).

To prepare 2% methylcellulose solution stock, for spheroid formation, 2 g methylcellulose is dissolved in 100 mL deionized water at 80°C. Then, the solution will be autoclaved then mixed overnight at 4°C. Solution can be kept at 4°C.

To prepare poly-HEMA stock solution to coat plates, 1.3 g of poly-HEMA is dissolved in 33 mL 99% ethanol by mixing overnight at 37°C. The solution is then filtered through a 0.22 µm filter and can be stored at room temperature (20°C). To coat plates, 50 µL of poly-HEMA stock solution is added to each well of a 96-well plate. Plates will be left to dry overnight at room temperature in an open tissue culture BSC hood for ventilation, and should be partially-closed instead of fully-closed. After drying, poly-HEMA-coated plates can be stored at 4°C for up to one month.

Seed cells at 4 different densities (25k, 50k, 75k, and 150k cells/cm²) in serum-free DMEM medium containing 0.75% methylcellulose, with 6 wells for each density. FBS is not added to the growth medium because it will promote cell attachment to the vessel wall and thus prevent spheroid formation. The methylcellulose will prevent excessive aggregation of cells in spheroid culture and maintain even spheroid size. Spheroids should be split every 3-4 days, and can be dissociated into single cells using Accutase (see section 6.2.c) and can be reseeded into new spheroids. See Appendix 12.2 for detailed protocol.

6.2.c. Spheroid Cell Dissociation

This protocol is based on Honeder *et al*, 2021 [86] and Han *et al*, 2020 [16], and it accounts for dissociation of 6 wells for seeding densities 1-4 (see section 6.2b). Spheroids can be dissociated after 3-4 days, and dissociated cells can be used to form new spheroids to determine the percent viability (see section 7.2). Firstly, 24 1.5 mL Eppendorf tubes are prepared and labeled accordingly, with 6 tubes for each seeding density group. Seeded contents from each well are transferred into their corresponding Eppendorf tubes and diluted with 480 µL PBS. Spheroids are collected via centrifugation at 800g for 15 min, and resuspended in 150 µL Accutase. Spheroids are incubated at 37°C for 10 minutes, and 850 µL serum-free DMEM is added to bring the final volume to 1 mL. See Appendix 12.7 for full protocol.

6.2.d. Lentiviral infection protocol

This protocol was developed by Hess Lab [87]. First, one should count and plate the number of cells listed above for each infection in a well of a 24-well plate (see Appendix 12.9). After filtering with a 0.45µm filter, 1-2 mL of virus(es) are added to each well and left to incubate for 24 hours at 37°C. After incubating, viral media is aspirated from each well and replaced with 1 mL of proper media. 48-72 hours later, the cells should be ready for analysis by fluorescence or selection with antibiotics. See Appendix 12.9 for detailed protocol.

6.2.e. High-throughput CRISPRi screen protocol

High-throughput, pooled CRISPRi-dCas9 loss-of-function screening protocol is based on Mathiowetz *et al*, 2023 [88] and may be modified to better suit this project's needs. Constitutive Cas9-expressing cells are first generated through lentiviral transduction. A dose response analysis is then performed to determine the concentration of cytotoxic compounds. For antibiotic selection, this analysis can either be a drug-resistance screen or a drug-sensitivity screen. Next, the dynamic range of the reporter for fluorescence-based assays is determined (aim for a greater distance between high and low fluorescence intensity bins). The sgRNA library is prepared, either by amplification and packaging into lentivirus or purchased as pre-packaged in lentivirus. The sgRNA is introduced via lentiviral transduction (note: this step can be done either together or separately from the Cas-9 introduction step). sgRNAs present in each cell population are then identified by extracting genomic DNA (gDNA) from frozen cell pellets and amplified by

polymerase chain reaction (PCR). The pooled sgRNA library is sent for deep sequencing and readings are analyzed using casTLE. See Appendix 12.10 for detailed protocol.

6.3. Final Prototype

6.3.a. A549

A549 is an adherent cell line isolated from the lung tissue of a White, 58-year-old male with lung cancer [28]. This NSCLC adenocarcinoma has a doubling time of 22 hours [28]. A549 exhibits mutations for TP53, EGFR, PIK3CA, ALK, and PTEN [27]. Since TP53 is the most frequently mutated gene in human cancer, the presence of this mutation is crucial for biological relevance to half of all cancer cases [57]. A549 is very sensitive to bleomycin (Z-score average of -1.04), somewhat sensitive to cisplatin (-0.12), and very sensitive to oxaliplatin (-1.45) [27], [32], [53]. Although one vial of these cells from ATCC is \$555.00, Hess lab provided 2 vials for this project [28]. Due to their 50 μm cell diameter, confluency in a T-75 flask was defined as 5,000,000 cells/10 mL.

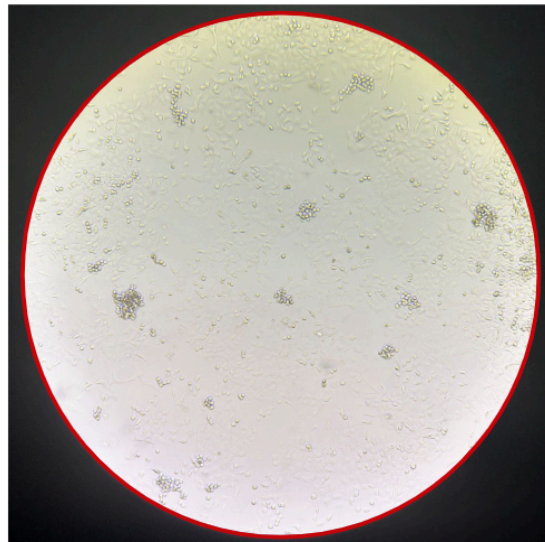


Figure 7: 10x Brightfield image of A549 Passage 5 cells taken at 200% confluency.

6.3.b. Poly-HEMA Plates

For full poly-HEMA stock solution plate preparation see Appendix 12.5.b. Briefly, add poly-HEMA to 99% ethanol, mix solution overnight at 37°C, filter solution into a glass bottle, and keep solution at room temperature until needed.

6.3.c. Spheroid Formation Protocol

For full materials and protocol see Appendix 12.5. Briefly, prepare methylcellulose stock solution, seed cells at multiple densities ranging from 200,000-150,000 cells/cm² (with 500 μL growth medium/cm²) in Poly-HEMA coated 96-well VWR plate, centrifuge spheroids, and leave in incubator for 3-4 days.

7. Testing

7.1. BioTek Cytation Imaging

In summary, Open BioTek Gen5 3.14 application, place plate in BioTek with A1 well in on the bottom right of the holder like pictured in Appendix 12.6, select new experiment and new

protocol making sure that the experiments and photos are saved in the same location and transferred together otherwise both will be lost, set protocol and imaging settings, and click run.

7.2. *Percent Viability after Spheroid Formation Determination*

The purpose of a spheroid cell dissociation experiment is to determine percent cell viability across seeding densities 1-4. Results from this experiment were meant to inform future percent live/dead experiments using mCherry for alive cells and Sytox Green for dead cells, with the aim to determine the seeding density showing a ~30% peak cell death rate during the initial growth phase of the spheroids [16]. Although this experiment was delayed to next semester due to an unforeseen maintenance issue with the Hess Lab's flow cytometry equipment, the protocol is outlined briefly here. After 1 mL of dissociated cells are obtained (see section 6.2.c), the cellular concentration and number of live and dead cells are counted using the CytoFLEX. The percent cell viability can be determined by comparing the number of live cells to the total number of cells (live and dead cells combined). See Appendix 12.8 for full protocol.

7.3. *ImageJ Analysis of BioTek Cytation Images*

To begin, choose the most in-focus image (figure 8a) for the desired well and upload it to an image analysis tool named ImageJ. Measure the scale bar in the corner of the image using the line tool and measurement feature. Set the scale under 'Analyze' → 'Set Scale' by setting the 'distance in pixels' to the measured amount of pixels from the images scale bar length. Then set the known length to the scale bar length in μm , and changing the 'Unit of length' to read μm in order to output results in micrometers. Select the 'Global' box to retain the scale bar incase of needing to restart the analysis process with an image, then select 'OK'.

Starting image analysis, change the image type to a 16-bit under 'Image' → 'Type' → '16-bit' in order to use the 'Threshold' feature. Next, select 'Image' → 'Adjust' → 'Threshold'. Uncheck the box titled 'Dark background' if selected, then move the ranges until the desired section shows the spheroids highlighted in red, as seen in figure 8b, then hit 'Apply'. Threshold recognizes pixels containing cells, or in this case spheroids. Next, go to 'Process' → 'Binary' → 'Watershed', which separates any spheroids that were close together but were recognized as one by the 'Threshold' feature. In order to block out any background noise from the edges of the well use the ellipse tool to encircle only the center of the well being analyzed (Figure 8c).

To analyze the highlighted spheroids click 'Analyze' → 'Analyze Particles'. Set the size range to 4147-infinity, as to ensure only spheroids of 20 or more cells were being analyzed. In the dropdown menu under 'Show:' select 'Outlines' then underneath select the 'Display results', 'Clear results', and 'Exclude on edges' boxes before hitting 'OK' (Figure 8d).

To save the results one of two windows will pop up. 1. A window asking to save measurements or 2. A results window. Should option 1 appear, select save measurements, but if option 2 appears, select 'File' → 'Save As' then save results as an excel file to a secure and known location. Combine the excel files into one master document for data analysis.

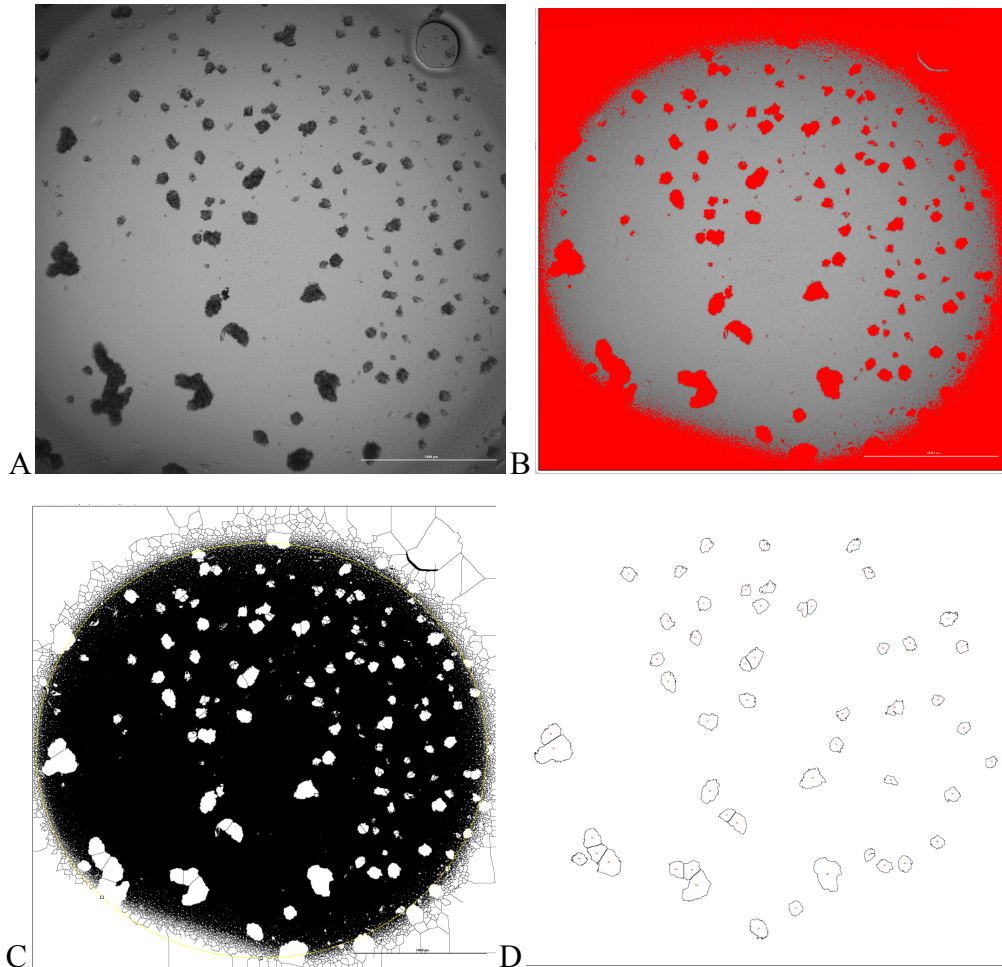


Figure 8: A. Initial image of well G1 from BioTek Cytation imager B. Threshold applied to image C. Binary→Watershed applied to the image with a thin yellow ellipse enclosing the well space D. Outlines of all spheroids analyzed for well G1

7.4. qPCR for Stemness Markers

This protocol was developed by Hess Lab [89]. A quantitative polymerase chain reaction (qPCR) will be performed to test for the expression level of stemness marker genes indicative of spheroid formation. For the purposes of this project, SOX2 and KLF4 will be the marker genes selected due to their association and upregulation in A549 spheroid cells [9].

First, RNA is extracted using the Qiagen RNeasy Kit. cDNA template is then synthesized using extracted RNA. Primers for gene(s) of interest are designed, along with preparing a loading control primer. qPCR is then performed using the Taqman method. CT values obtained for experimental samples compared to control are analyzed to determine if the expression of gene(s) of interest has been perturbed. See Appendix 12.11 for detailed protocol.

8. Results

The results from the ImageJ analysis were compiled into one master document to further analyze. The average number of spheroids per density and the average spheroid size per density were taken.

First the average spheroid size per density was taken by averaging all the data from the six wells of each density. There appeared to be no apparent correlation between size and density as the average was 5675, 8558, 6301, and 9232 μm^2 for 25k, 50k, 75k, and 150k cells/cm² respectively.

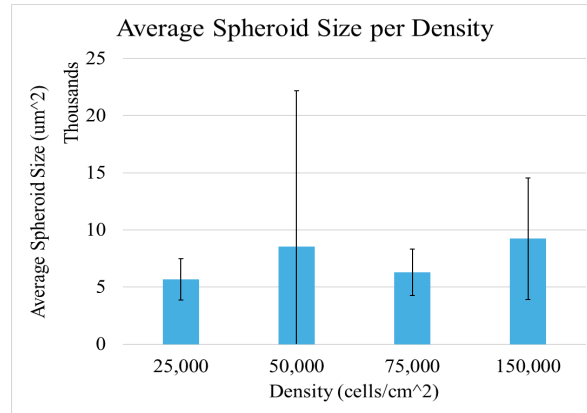


Figure 9: Average Spheroid Size per Density

Secondly, the average number of spheroids per density was taken by once again averaging all the data from the six wells of each density. There appeared to be a positive correlation between an increase of spheroid number to an increase in density. The average were 4, 23, 39, and 60 spheroids for 25k, 50k, 75k, and 150k cells/cm² respectively.

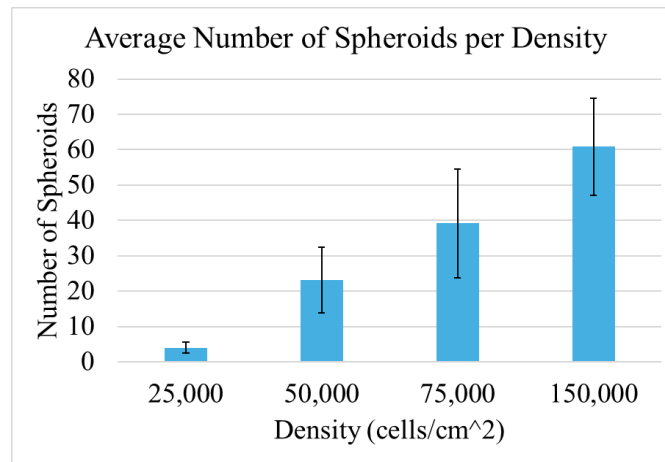


Figure 10: Average Number of Spheroids per Density.

9. Discussion

The passaging protocol for A549 cells was optimized to reflect the established doubling time of approximately 23 hours. This optimization ensured consistent cell viability and growth, avoiding over-confluence that could disrupt spheroid formation and cellular responses. By implementing a schedule of one passage every 48-72 hours, the team maintained cellular integrity, aligning with literature values and reducing variability caused by inconsistent handling. During the spheroid formation process, the results demonstrated a linear increase in spheroid number with seeding density. The average spheroid size remained within the 50–800 μm diameter range reported in the literature. Towards the end of the semester, the team could have made more spheroids to further back the results seen in figure 10, but issues with the CytoFLEX led to a change in plans. Next semester a repeat trial of seeding density variations using a

live/dead cytometry assay will confirm the preliminary findings regarding spheroid viability and optimal seeding conditions. Additionally, in creating uniform spheroids, the methylcellulose concentration for each well is then considered optimized.

While the team was successful in creating spheroids, performing qPCR for SOX2 will have to be done next semester. It is hypothesized that gene expression differences are to be seen when comparing 3D spheroids and 2D monolayer cultures, with the idea that a 3D model is better at mimicking *in vivo* tumor microenvironments. SOX2 amplifies the ability of cancer cells to sustain growth and resist apoptosis, which is why it is a significant measurement of tumorigenicity. The SOX2 qPCR would provide valuable information pre-CRISPRi screen. For instance, the CRISPRi screen could knockdown a gene such as SOX2, allowing for repression at the transcriptional level to see how the gene played a role in growth and survival. The more favorable a gene is, the more DNA damage will be seen due to its repression. While in hindsight this process seems to be moving in the right direction, there is currently no data on CRISPRi screening within a A549 spheroid model, which leaves the results and processes involved open to discussion.

Lastly, the final future work aspect of this project involves testing the processes of dissociation, staining, and fixing the spheroid cells for γ H2AX. This will help determine how the specific genes targeted respond to the repression from a CRISPRi screen. γ H2AX is a marker for DSBs, which will either be increased or decreased as a result of gene repression. DSBs are one of the more severe forms of DNA damage, where both DNA helix strands are broken simultaneously. These breaks can be the result in many cases from ionizing radiation, reactive oxygen species, or chemotherapeutic agents for cancer treatment. By utilizing γ H2AX staining in a 3D environment, more accurate predictions can be made about what genes contribute to tumor progression and metastasis. The goal is to refine these methods and achieve these significant milestones by May 2025.

10. Conclusions

This semester the team was able to select a viable adenocarcinoma non-small cell lung cancer cell line: A549. The cell line was chosen based on characteristics such as: doubling time, adherence properties, and key tumor protein 53 mutation. The cell line was maintained and counted to stay at a consistent 60-80% confluency, which for A549 cells is 3 to 4 million cells per vial. For spheroid formation, poly-HEMA treated plates were chosen due to their cost-effectiveness and ability to produce a consistent size and shape of spheroids. The crosslinking of poly-HEMA prevents cell attachment and promotes cell-cell aggregation into 3D spheroids that are better at mimicking an *in vivo* environment than a traditional 2D structure. Current research has an optimal spheroid size for the A549 cell line around 500 μ m, so testing data from next semester will be used to compare against this value. By optimizing a spheroid formation protocol, the team will be able to begin working on a γ H2AX staining protocol at the beginning of next semester.

1. References

- [1] "Cancer Facts & Figures 2024." Accessed: Dec. 12, 2024. [Online]. Available: <https://www.cancer.org/research/cancer-facts-statistics/all-cancer-facts-figures/2024-cancer-facts-figures.html>
- [2] M. Kapałczyńska *et al.*, "2D and 3D cell cultures – a comparison of different types of

- cancer cell cultures,” *Arch. Med. Sci. AMS*, vol. 14, no. 4, pp. 910–919, Jun. 2018, doi: 10.5114/aoms.2016.63743.
- [3] “2D vs 3D cell culture | Learning Center,” 2D vs 3D cell culture | Learning Center | UPM Biomedicals. Accessed: Oct. 03, 2024. [Online]. Available: <https://www.upmbiomedicals.com/resource-center/learning-center/what-is-3d-cell-culture/2d-versus-3d-cell-culture/>
- [4] K. Juárez-Moreno, D. Chávez-García, G. Hirata, and R. Vazquez-Duhalt, “Monolayer (2D) or spheroids (3D) cell cultures for nanotoxicological studies? Comparison of cytotoxicity and cell internalization of nanoparticles,” *Toxicol. In Vitro*, vol. 85, p. 105461, Dec. 2022, doi: 10.1016/j.tiv.2022.105461.
- [5] “2D vs 3D cell culture | Learning Center | UPM Biomedicals.” Accessed: Oct. 13, 2024. [Online]. Available: <https://www.upmbiomedicals.com/resource-center/learning-center/what-is-3d-cell-culture/2d-versus-3d-cell-culture/>
- [6] admin.facellitate, “What are primary cells : advantages and limitations,” faCellitate. Accessed: Oct. 14, 2024. [Online]. Available: <https://facellitate.com/what-are-primary-cells-advantages-and-limitations/>
- [7] J.-P. Gillet, S. Varma, and M. M. Gottesman, “The Clinical Relevance of Cancer Cell Lines,” *JNCI J. Natl. Cancer Inst.*, vol. 105, no. 7, pp. 452–458, Apr. 2013, doi: 10.1093/jnci/djt007.
- [8] C. He *et al.*, “CRISPR screen in cancer: status quo and future perspectives,” *Am. J. Cancer Res.*, vol. 11, no. 4, pp. 1031–1050, Apr. 2021.
- [9] C. Choe, H. Kim, S. Min, S. Park, J. Seo, and S. Roh, “SOX2, a stemness gene, induces progression of NSCLC A549 cells toward anchorage-independent growth and chemoresistance to vinblastine,” *OncoTargets Ther.*, vol. 11, pp. 6197–6207, Sep. 2018, doi: 10.2147/OTT.S175810.
- [10] S. Dasari and P. B. Tchounwou, “Cisplatin in cancer therapy: molecular mechanisms of action,” *Eur. J. Pharmacol.*, vol. 740, pp. 364–378, Oct. 2014, doi: 10.1016/j.ejphar.2014.07.025.
- [11] B. Devanabanda and A. Kasi, “Oxaliplatin,” in *StatPearls*, Treasure Island (FL): StatPearls Publishing, 2024. Accessed: Oct. 14, 2024. [Online]. Available: <http://www.ncbi.nlm.nih.gov/books/NBK557690/>
- [12] “Bleomycin - StatPearls - NCBI Bookshelf.” Accessed: Oct. 14, 2024. [Online]. Available: <https://www.ncbi.nlm.nih.gov/books/NBK555895/>
- [13] L. Morrison, J. A. Laukkanen, K. Ronkainen, S. Kurl, J. Kauhanen, and A. T. Toriola, “Inflammatory biomarker score and cancer: A population-based prospective cohort study,” *BMC Cancer*, vol. 16, p. 80, Feb. 2016, doi: 10.1186/s12885-016-2115-6.
- [14] J. Fan *et al.*, “Multiplex gene quantification as digital markers for extremely rapid evaluation of chemo-drug sensitivity,” *Patterns*, vol. 2, no. 10, p. 100360, Sep. 2021, doi: 10.1016/j.patter.2021.100360.
- [15] N. Takahashi *et al.*, “3D Culture Models with CRISPR Screens Reveal Hyperactive NRF2 as a Prerequisite for Spheroid Formation via Regulation of Proliferation and Ferroptosis,” *Mol. Cell*, vol. 80, no. 5, pp. 828–844.e6, Dec. 2020, doi: 10.1016/j.molcel.2020.10.010.
- [16] K. Han *et al.*, “CRISPR screens in cancer spheroids identify 3D growth-specific vulnerabilities,” *Nature*, vol. 580, no. 7801, pp. 136–141, Apr. 2020, doi: 10.1038/s41586-020-2099-x.

- [17] J. R. Ferrarone *et al.*, “Genome-wide CRISPR screens in spheroid culture reveal that the tumor suppressor LKB1 inhibits growth via the PIKFYVE lipid kinase,” *Proc. Natl. Acad. Sci.*, vol. 121, no. 21, p. e2403685121, May 2024, doi: 10.1073/pnas.2403685121.
- [18] T. Stiff, S. Bayraktar, P. Dama, J. Stebbing, and L. Castellano, “CRISPR screens in 3D tumourspheres identified miR-4787-3p as a transcriptional start site miRNA essential for breast tumour-initiating cell growth,” *Commun. Biol.*, vol. 7, no. 1, pp. 1–11, Jul. 2024, doi: 10.1038/s42003-024-06555-1.
- [19] S. Gopal, A. L. Rodrigues, and J. S. Dordick, “Exploiting CRISPR Cas9 in Three-Dimensional Stem Cell Cultures to Model Disease,” *Front. Bioeng. Biotechnol.*, vol. 8, Jun. 2020, doi: 10.3389/fbioe.2020.00692.
- [20] “Hess, Gaelen,” Integrated Program in Biochemistry. Accessed: Oct. 05, 2024. [Online]. Available: <https://ipib.wisc.edu/staff/hess-gaelen/>
- [21] R. Blanco *et al.*, “A Gene-Alteration Profile of Human Lung Cancer Cell Lines,” *Hum. Mutat.*, vol. 30, no. 8, pp. 1199–1206, Aug. 2009, doi: 10.1002/humu.21028.
- [22] CDC, “Biosafety in Microbiological and Biomedical Laboratories—6th Edition,” p. 37, Jun. 2020.
- [23] “Freezing Cells - US.” Accessed: Oct. 13, 2024. [Online]. Available: <https://www.thermofisher.com/us/en/home/references/gibco-cell-culture-basics/cell-culture-protocols/freezing-cells.html>
- [24] M. Selby *et al.*, “3D Models of the NCI60 Cell Lines for Screening Oncology Compounds,” *SLAS Discov. Adv. Life Sci. R D*, vol. 22, no. 5, pp. 473–483, Jun. 2017, doi: 10.1177/2472555217697434.
- [25] “NCI-H23 [H23] - CRL-5800 | ATCC.” Accessed: Sep. 25, 2024. [Online]. Available: <https://www.atcc.org/products/crl-5800>
- [26] N. Takahashi *et al.*, “3D Culture Models with CRISPR Screens Reveal Hyperactive NRF2 as a Prerequisite for Spheroid Formation via Regulation of Proliferation and Ferroptosis,” *Mol. Cell*, vol. 80, no. 5, pp. 828-844.e6, Dec. 2020, doi: 10.1016/j.molcel.2020.10.010.
- [27] “Cell Line:905942 - Cancerrxgene - Genomics of Drug Sensitivity in Cancer.” Accessed: Sep. 27, 2024. [Online]. Available: <https://www.cancerrxgene.org/cellline/NCI-H23/905942>
- [28] “A549 - CCL-185 | ATCC.” Accessed: Sep. 25, 2024. [Online]. Available: <https://www.atcc.org/products/ccl-185>
- [29] L. Korrodi-Gregório, V. Soto-Cerrato, R. Vitorino, M. Fardilha, and R. Pérez-Tomás, “From Proteomic Analysis to Potential Therapeutic Targets: Functional Profile of Two Lung Cancer Cell Lines, A549 and SW900, Widely Studied in Pre-Clinical Research,” *PLoS ONE*, vol. 11, no. 11, p. e0165973, Nov. 2016, doi: 10.1371/journal.pone.0165973.
- [30] “OVCAR-5 Human Cancer Cell Line | SCC259.” Accessed: Sep. 27, 2024. [Online]. Available: https://www.emdmillipore.com/US/en/product/OVCAR-5-Human-Cancer-Cell-Line,MM_NF-SCC259?ReferrerURL=https%3A%2F%2Fwww.google.com%2F
- [31] B. C. Baguley *et al.*, “Resistance mechanisms determining the in vitro sensitivity to paclitaxel of tumour cells cultured from patients with ovarian cancer,” *Eur. J. Cancer Oxf. Engl.* 1990, vol. 31A, no. 2, pp. 230–237, 1995, doi: 10.1016/0959-8049(94)00472-h.
- [32] “Cell Line:905969 - Cancerrxgene - Genomics of Drug Sensitivity in Cancer.” Accessed: Sep. 27, 2024. [Online]. Available: <https://www.cancerrxgene.org/cellline/OVCAR-5/905969>
- [33] “OVCAR-5/RFP cell line.” Accessed: Sep. 30, 2024. [Online]. Available:

- <https://www.mybiosource.com/cell-line/ovcar-5-rfp/168983>
- [34] W. Asghar *et al.*, “In Vitro Three-Dimensional Cancer Culture Models,” *Cancer Target. Drug Deliv. Elus. Dream*, pp. 635–665, Jul. 2013, doi: 10.1007/978-1-4614-7876-8-24.
- [35] L. Carroll, B. Tiwari, J. Curtin, and J. Wanigasekara, “U-251MG Spheroid Generation Using Hanging Drop Method Protocol,” May 2021, Accessed: Sep. 27, 2024. [Online]. Available: <https://www.protocols.io/view/u-251mg-spheroid-generation-using-hanging-drop-met-btstn-nen>
- [36] “Generation and Cultivation of 3D Cell Models,” *ibidi*. Accessed: Oct. 10, 2024. [Online]. Available: <https://ibidi.com/content/1037-generation-and-cultivation-of-3d-cell-models>
- [37] A. P. P. Guimaraes, I. R. Calori, H. Bi, and A. C. Tedesco, “SpheroMold: modernizing the hanging drop method for spheroid culture,” *Front. Drug Deliv.*, vol. 4, Jun. 2024, doi: 10.3389/fddev.2024.1397153.
- [38] “Perfecta3D Hanging Drop Plates - 96-Well.” Accessed: Oct. 10, 2024. [Online]. Available: <https://www.selectscience.net/product/perfecta3d-hanging-drop-plates-96-well>
- [39] S. S. Nazari, “Generation of 3D Tumor Spheroids with Encapsulating Basement Membranes for Invasion Studies,” *Curr. Protoc. Cell Biol.*, vol. 87, no. 1, p. e105, 2020, doi: 10.1002/cpcb.105.
- [40] M. A. Badea *et al.*, “Influence of Matrigel on Single- and Multiple-Spheroid Cultures in Breast Cancer Research,” *SLAS Discov.*, vol. 24, no. 5, pp. 563–578, Jun. 2019, doi: 10.1177/2472555219834698.
- [41] “Overview of Cell Viability and Survival,” Cell Signaling Technology. Accessed: Sep. 27, 2024. [Online]. Available: <https://www.cellsignal.com/science-resources/cell-viability-and-survival>
- [42] M. M. Silva, C. R. R. Rocha, G. S. Kinker, A. L. Pelegrini, and C. F. M. Menck, “The balance between NRF2/GSH antioxidant mediated pathway and DNA repair modulates cisplatin resistance in lung cancer cells,” *Sci. Rep.*, vol. 9, p. 17639, Nov. 2019, doi: 10.1038/s41598-019-54065-6.
- [43] A. R. Hyler *et al.*, “Fluid shear stress impacts ovarian cancer cell viability, subcellular organization, and promotes genomic instability,” *PLoS ONE*, vol. 13, no. 3, p. e0194170, Mar. 2018, doi: 10.1371/journal.pone.0194170.
- [44] R. Zhang, X. Chen, S. Fu, L. Xu, and J. Lin, “A small molecule STAT3 inhibitor, LLL12, enhances cisplatin- and paclitaxel-mediated inhibition of cell growth and migration in human ovarian cancer cells,” *Oncol. Rep.*, vol. 44, no. 3, pp. 1224–1232, Sep. 2020, doi: 10.3892/or.2020.7667.
- [45] J. N. Søndergaard, K. Geng, C. Sommerauer, I. Atanasoai, X. Yin, and C. Kutter, “Successful delivery of large-size CRISPR/Cas9 vectors in hard-to-transfect human cells using small plasmids,” *Commun. Biol.*, vol. 3, p. 319, Jun. 2020, doi: 10.1038/s42003-020-1045-7.
- [46] T. Kim, I. Doh, and Y.-H. Cho, “On-chip three-dimensional tumor spheroid formation and pump-less perfusion culture using gravity-driven cell aggregation and balanced droplet dispensing,” *Biomicrofluidics*, vol. 6, no. 3, p. 034107, Jul. 2012, doi: 10.1063/1.4739460.
- [47] C. D. House, L. Hernandez, and C. M. Annunziata, “In vitro Enrichment of Ovarian Cancer Tumor-initiating Cells,” *J. Vis. Exp. JoVE*, no. 96, p. 52446, Feb. 2015, doi: 10.3791/52446.
- [48] “OVCAR-5 Cell Avalanche® Transfection Reagent – EZ Biosystems.” Accessed: Sep. 25, 2024. [Online]. Available:

- <https://ezbiosystems.com/product/ovcar-5-cell-avalanche-transfection-reagent/>
- [49] J.-C. Cheng, C. Klausen, and P. C. K. Leung, “Hypoxia-inducible factor 1 alpha mediates epidermal growth factor-induced down-regulation of E-cadherin expression and cell invasion in human ovarian cancer cells,” *Cancer Lett.*, vol. 329, no. 2, pp. 197–206, Feb. 2013, doi: 10.1016/j.canlet.2012.10.029.
- [50] “Cell Doubling Time Calculator.” Accessed: Sep. 25, 2024. [Online]. Available: <https://www.omnicalculator.com/biology/cell-doubling-time>
- [51] “Cellosaurus - Cell line encyclopedia.” Accessed: Sep. 25, 2024. [Online]. Available: <https://www.cellosaurus.org/index.html>
- [52] “Home page - Cancerrxgene - Genomics of Drug Sensitivity in Cancer.” Accessed: Sep. 27, 2024. [Online]. Available: <https://www.cancerrxgene.org/>
- [53] “Cell Line:905949 - Cancerrxgene - Genomics of Drug Sensitivity in Cancer.” Accessed: Sep. 27, 2024. [Online]. Available: <https://www.cancerrxgene.org/cellline/A549/905949>
- [54] “TP53 tumor protein p53 BCC7 LFS1 P53 TRP53 | Sigma-Aldrich.” Accessed: Sep. 19, 2024. [Online]. Available: <https://www.sigmaaldrich.com/US/en/genes/tp53>
- [55] “EGFR epidermal growth factor receptor ERBB ERBB1 HER1 PIG61 mENA | Sigma-Aldrich.” Accessed: Sep. 20, 2024. [Online]. Available: <https://www.sigmaaldrich.com/US/en/genes/egfr>
- [56] “KRAS Kirsten rat sarcoma viral oncogene homolog C-K-RAS CFC2 K-RAS2A K-RAS2B K-RAS4A K-RAS4B KI-RAS KRAS1 KRAS2 NS NS3 RASK2 | Sigma-Aldrich.” Accessed: Sep. 19, 2024. [Online]. Available: <https://www.sigmaaldrich.com/US/en/genes/kras>
- [57] R. Mirgayazova *et al.*, “Therapeutic Editing of the TP53 Gene: Is CRISPR/Cas9 an Option?,” *Genes*, vol. 11, no. 6, p. 704, Jun. 2020, doi: 10.3390/genes11060704.
- [58] “PIK3CA phosphatidylinositol-4,5-bisphosphate 3-kinase, catalytic subunit alpha CLOVE CWS5 MCAP MCM MCMT C PI3K p110-alpha | Sigma-Aldrich.” Accessed: Sep. 27, 2024. [Online]. Available: <https://www.sigmaaldrich.com/US/en/genes/pik3ca>
- [59] S.-Q. Hou, M. Ouyang, A. Brandmaier, H. Hao, and W. H. Shen, “PTEN in the maintenance of genome integrity: From DNA replication to chromosome segregation,” *BioEssays News Rev. Mol. Cell. Dev. Biol.*, vol. 39, no. 10, p. 10.1002/bies.201700082, Oct. 2017, doi: 10.1002/bies.201700082.
- [60] E. Ó hAinmhire, S. M. Quartuccio, W. Cheng, R. A. Ahmed, S. M. King, and J. E. Burdette, “Mutation or Loss of p53 Differentially Modifies TGFβ Action in Ovarian Cancer,” *PLoS ONE*, vol. 9, no. 2, p. e89553, Feb. 2014, doi: 10.1371/journal.pone.0089553.
- [61] “CLDN4 claudin 4 [Homo sapiens (human)] - Gene - NCBI.” Accessed: Sep. 27, 2024. [Online]. Available: <https://www.ncbi.nlm.nih.gov/gene/1364>
- [62] T. Riss and O. J. Trask, “Factors to consider when interrogating 3D culture models with plate readers or automated microscopes,” *In Vitro Cell. Dev. Biol. Anim.*, vol. 57, no. 2, pp. 238–256, 2021, doi: 10.1007/s11626-020-00537-3.
- [63] H. Shi *et al.*, “3-Dimensional mesothelioma spheroids provide closer to natural pathophysiological tumor microenvironment for drug response studies,” *Front. Oncol.*, vol. 12, p. 973576, Aug. 2022, doi: 10.3389/fonc.2022.973576.
- [64] Y.-C. Chen, X. Lou, Z. Zhang, P. Ingram, and E. Yoon, “High-Throughput Cancer Cell Sphere Formation for Characterizing the Efficacy of Photo Dynamic Therapy in 3D Cell Cultures,” *Sci. Rep.*, vol. 5, no. 1, p. 12175, Jul. 2015, doi: 10.1038/srep12175.

- [65] A. P. Aijian and R. L. Garrell, “Digital Microfluidics for Automated Hanging Drop Cell Spheroid Culture,” *SLAS Technol.*, vol. 20, no. 3, pp. 283–295, Jun. 2015, doi: 10.1177/2211068214562002.
- [66] R. Teixeira Polez, N. Huynh, C. S. Pridgeon, J. J. Valle-Delgado, R. Harjumäki, and M. Österberg, “Insights into spheroids formation in cellulose nanofibrils and Matrigel hydrogels using AFM-based techniques,” *Mater. Today Bio*, vol. 26, p. 101065, Jun. 2024, doi: 10.1016/j.mtbio.2024.101065.
- [67] N.-E. Ryu, S.-H. Lee, and H. Park, “Spheroid Culture System Methods and Applications for Mesenchymal Stem Cells,” *Cells*, vol. 8, no. 12, p. 1620, Dec. 2019, doi: 10.3390/cells8121620.
- [68] “Low Attachment Cell Culture Plates, Dishes, and Flasks - US.” Accessed: Sep. 27, 2024. [Online]. Available: <https://www.thermofisher.com/us/en/home/life-science/cell-culture/organoids-spheroids-3d-cell-culture/3d-cell-culture-plates-dishes-flasks.html>
- [69] “Trypsin-EDTA (0.05%), phenol red.” Accessed: Sep. 27, 2024. [Online]. Available: <https://www.thermofisher.com/order/catalog/product/25300054>
- [70] “RNase A, DNase and protease-free (10 mg/mL).” Accessed: Sep. 27, 2024. [Online]. Available: <https://www.thermofisher.com/order/catalog/product/EN0531>
- [71] R. Foty, “A Simple Hanging Drop Cell Culture Protocol for Generation of 3D Spheroids,” *J. Vis. Exp. JoVE*, no. 51, p. 2720, May 2011, doi: 10.3791/2720.
- [72] “Corning® Matrigel® Basement Membrane Matrix -.” Accessed: Sep. 27, 2024. [Online]. Available: <https://www.sigmaaldrich.com/US/en/substance/corningmatrigelbasementmembranematrix1234598765>
- [73] “A549 Cell Line Maintenance · Benchling.” Accessed: Oct. 05, 2024. [Online]. Available: https://benchling.com/uw_hesslab/f/lib_tD5E8G7L-cell-line-maintenance/etr_5Skd74TL-a549-cell-line-maintenance/edit
- [74] “Growth rate and doubling time – Cellseeker Inventory.” Accessed: Dec. 12, 2024. [Online]. Available: <https://www.cellseeker.org/cellcalc/growth-rate-and-doubling-time/>
- [75] “T75 Flasks - US.” Accessed: Oct. 13, 2024. [Online]. Available: <https://www.thermofisher.com/us/en/home/life-science/cell-culture/cell-culture-plastics/cell-culture-flasks/t75-flasks.html>
- [76] “DMEM, high glucose.” Accessed: Dec. 17, 2024. [Online]. Available: <https://www.thermofisher.com/order/catalog/product/11965118>
- [77] “Penicillin-Streptomycin (5,000 U/mL).” Accessed: Oct. 13, 2024. [Online]. Available: <https://www.thermofisher.com/order/catalog/product/15070063>
- [78] “Fetal Bovine Serum, Value (formerly USDA-approved in North America or qualified, Brazil in other regions).” Accessed: Oct. 13, 2024. [Online]. Available: <https://www.thermofisher.com/order/catalog/product/A5256701>
- [79] “PBS, pH 7.4.” Accessed: Oct. 13, 2024. [Online]. Available: <https://www.thermofisher.com/order/catalog/product/10010023>
- [80] “VWR® Polystyrene Assay Plates, 96- and 384-Well,” VWR. Accessed: Dec. 12, 2024. [Online]. Available: <https://us.vwr.com/store/product/42378328/vwr-polystyrene-assay-plates-96-and-384-well>
- [81] “Ethanol, Absolute (200 Proof), Molecular Biology Grade, Fisher BioReagents, Quantity: 100 mL | Fisher Scientific.” Accessed: Dec. 12, 2024. [Online]. Available:

- <https://www.fishersci.com/shop/products/ethanol-absolute-200-proof-molecular-biology-grade-fisher-bioreagents-3/BP28184>
- [82] “Poly(2-hydroxyethyl methacrylate) BioReagent, powder, cell culture mammalian 25249-16-5.” Accessed: Dec. 12, 2024. [Online]. Available: <https://www.sigmaaldrich.com/US/en/product/sigma/p3932>
- [83] “Methyl cellulose viscosity: 4,000 cP 9004-67-5.” Accessed: Dec. 12, 2024. [Online]. Available: <https://www.sigmaaldrich.com/US/en/product/sigma/m0512>
- [84] “StemPro™ Accutase™ Cell Dissociation Reagent.” Accessed: Dec. 12, 2024. [Online]. Available: <https://www.thermofisher.com/order/catalog/product/A1110501>
- [85] “Passaging Cells · Benchling.” Accessed: Oct. 05, 2024. [Online]. Available: https://benchling.com/uw_hesslab/f/lib_jUyOwxBE-tissue-culture-protocols/prt_PWYXyJOG-passaging-cells/edit
- [86] S. Honeder *et al.*, “Adipose Triglyceride Lipase Loss Promotes a Metabolic Switch in A549 Non–Small Cell Lung Cancer Cell Spheroids,” *Mol. Cell. Proteomics MCP*, vol. 20, p. 100095, May 2021, doi: 10.1016/j.mcpro.2021.100095.
- [87] “Small scale infection of mammalian cells with Lentivirus: Adherent Cells Version · Benchling.” Accessed: Oct. 05, 2024. [Online]. Available: https://benchling.com/uw_hesslab/f/lib_jUyOwxBE-tissue-culture-protocols/prt_Sdt7mJ5N-small-scale-infection-of-mammalian-cells-with-lentivirus-adherent-cells-version/edit
- [88] A. J. Mathiowetz, M. A. Roberts, D. W. Morgens, J. A. Olzmann, and Z. Li, “Protocol for performing pooled CRISPR-Cas9 loss-of-function screens,” *STAR Protoc.*, vol. 4, no. 2, p. 102201, Mar. 2023, doi: 10.1016/j.xpro.2023.102201.
- [89] “qPCR Protocol · Benchling.” Accessed: Oct. 10, 2024. [Online]. Available: https://benchling.com/uw_hesslab/f/lib_xtPaCTIU-extraction-and-staining-protocols/prt_Va3aG4Cc-qpcr-protocol/edit